

Testicular Implants

Silicone Gel-Filled
Highly Polymerized Silicone

POLYTECH

Health & Aesthetics



Silicone gel-filled and highly polymerized silicone testicular implants are available in up to 5 sizes.

Testicular Implants, € 0483

Product Description

Testicular implants with smooth surface are available in two different types, which are both suitable for an attachment to the tissue.

The **silicone-gel filled** implants are fitted with a short piece of removable teflon to guide the needle and therefore to prevent a puncture of the implant.

The **highly polymerized** testicular implants have a reinforced cap that allows the attachment of the implant with a suture through the cap.

Information on Materials

The shell, which contains a barrier layer, is filled with either a highly cross-linked silicone gel or a highly polymerized silicone. Both types of testicular implants have a smooth surface.

Indications

- Reconstruction after orchidectomy
- Reconstruction after trauma
- Correction of asymmetry
- Congenital defects
- Implant replacement

Important Information

Testicular implants are delivered individually packed and sterile.

The teflon needle guide of the silicone gel-filled implants must always be removed before wound closure!

Attention: Please note the detailed instructions for use included with each product.




Manufacturer:

POLYTECH Health & Aesthetics GmbH
 Altheimer Str. 32 • 64807 Dieburg • Germany
 Tel.: +49 (0) 60 71 98 63 0 • Fax: +49 (0) 60 71 98 63 30
 www.polytech-health-aesthetics.com • E-mail: info@polytechhealth.com

Order no.	Weight (g)	Dimensions (mm)	
		A	B


Testicular Implants, silicone gel-filled, smooth

10100-005	5	20	22
10100-012	12	26	33
10100-021	21	32	42
10100-027	27	33	48
10100-036	36	39	50

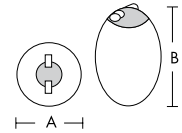

Testicular Implants, highly polymerized silicone, smooth

10101-005	5	20	22
10101-012	12	26	33
10101-021	21	32	42
10101-027	27	33	48
10101-036	36	39	50

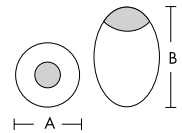

Approximate measurements