

A

Blades
for Large-Bone Power Systems

A-2

Blades
for Small-Bone Power Systems

A-15

Sternum Saw Blades

A-31

Burs

A-34

Rasps

A-40

K-wires + Steinmann Pins

A-41

Twist Drills + Taps

A-43

MICROAIRE® and Linvatec®

MICROAIRE®
For Surgery. For Life.™

Blades, Burs, Rasps, K-wires, Steinmann Pins, Twist Drills and Taps

BLADES for Large-Bone Power Systems

MicroAire® was the first company to offer double-sterile-packaged replacement large-power blades, and our commitment to value pricing and first-class service continues today.

For use with the following instruments:

MicroAire® Series 7000

Hall® PowerPro®, Versipower® and Series 4, PowerPro, M-Power, M-Power²



StableCut®

- Patented tooth design
- Aggressive angles cut along a straight path
- Increased ejection of debris
- Decreased friction reduces bone necrosis and implant failure
- Top choice of implant companies and orthopedic surgeons



Aggressive Offset

- Tooth geometry similar to PhotoEtch Aggressive blade
- Two rows of teeth formed by bending alternating teeth away from the surface of the blade



PhotoEtch Aggressive

- Stainless steel photo-chemically etched into a dual-thickness blade
- Two rows of aggressive-angled teeth
- Dual thickness reduces friction when used with a slotted cutting guide



Standard Offset

- Tooth geometry similar to Diamond-Ground blades
- Two rows of teeth formed by bending alternating teeth away from the surface of the blade



Diamond Ground

- Tooth pattern is cut into blade
- Angles are ground into the faces of teeth
- Two rows of sharp teeth
- Cutting thickness is same as the blade itself



D-Cut Replacement

- Small "double-teeth" on a straight path instead of a radius
- Channels between teeth help to move material away from cutting surface



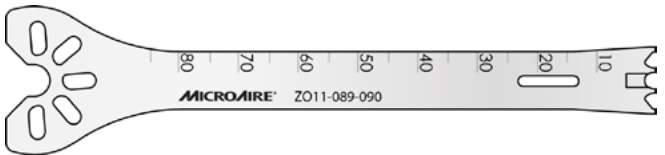
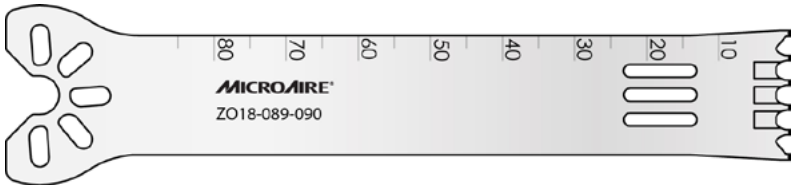
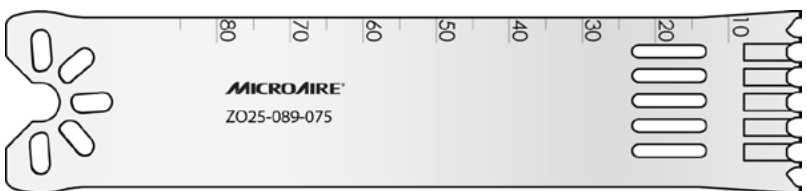
Reciprocating Blades

- Used in total hip and total joint arthroplasty
- Available in both single- and double-edged patterns

Oscillating Saw Blades

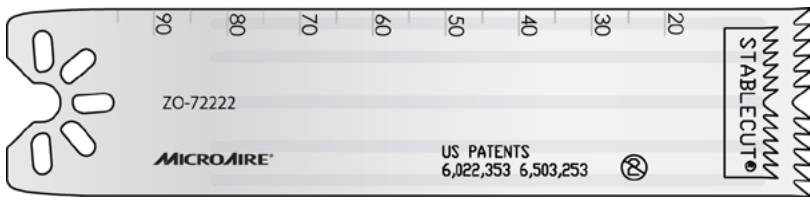
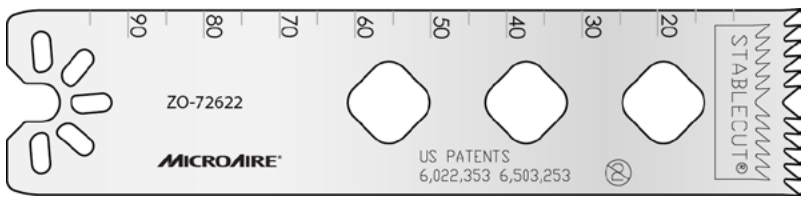
MicroAire® and Hall® Large-Power Saws

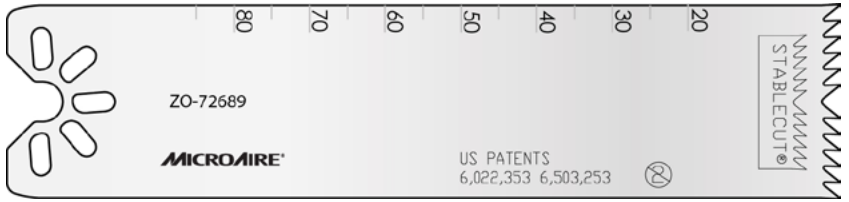
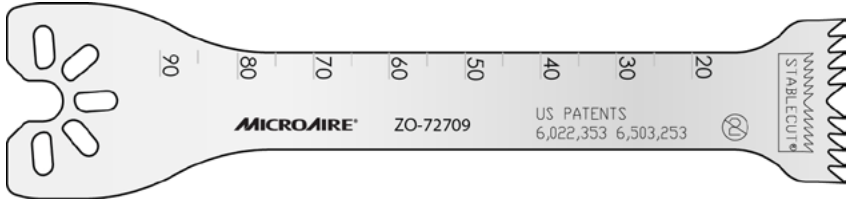
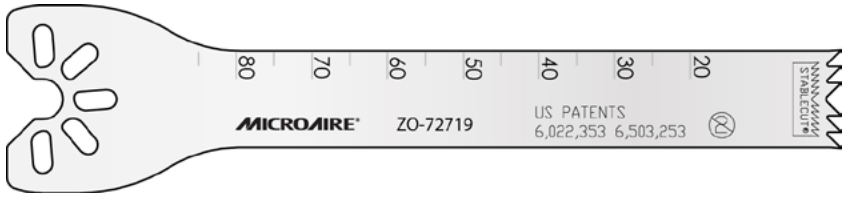
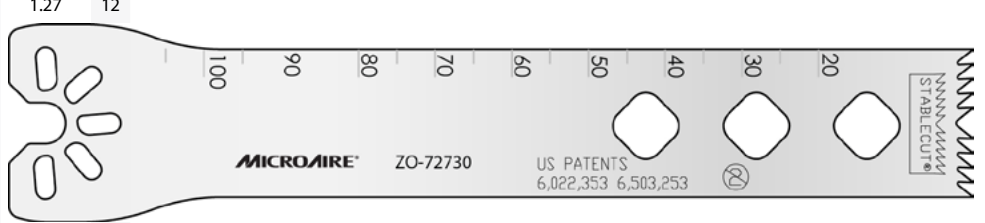
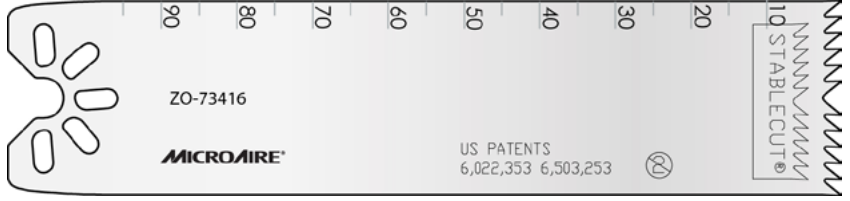
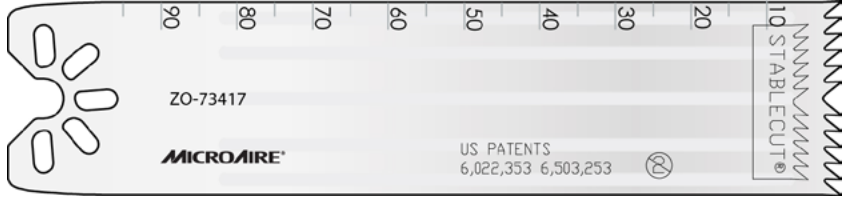
D-Cut®

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
Z011-089-090	11	90	0.89	0.89	16	
Z011-119-090	11	90	1.19	1.19	16	
Z011-127-090	11	90	1.27	1.27	16	
Z018-089-090	18	90	0.89	0.89	16	
Z018-127-090	18	90	1.27	1.27	16	
Z018-127-100	18	100	1.27	1.27	16	
Z018-135-090	18	90	1.35	1.35	16	
Z025-089-075	25	75	0.89	0.89	16	
Z025-089-090	25	90	0.89	0.89	16	
Z025-107-090	25	90	1.07	1.07	16	
Z025-119-090	25	90	1.19	1.19	16	
Z025-127-090	25	90	1.27	1.27	16	
Z025-127-100	25	100	1.27	1.27	16	
Z025-135-090	25	90	1.35	1.35	16	
Z025-147-090	25	90	1.47	1.47	16	

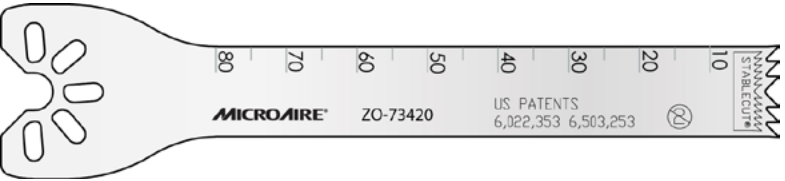
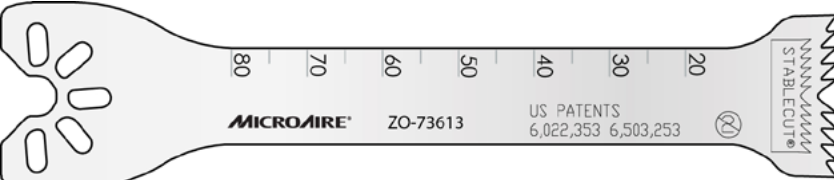
MicroAire® and Hall® Large-Power Saws

StableCut®

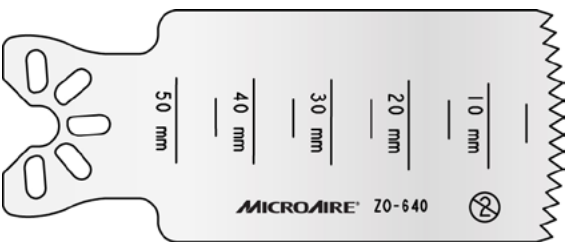
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
Z0-72222	25	90	1.27	1.27	12	
Z0-72622	25	90	1.19	1.19	12	
Z0-72623	25	90	1.27	1.27	12	
Z0-72625	25	90	1.37	1.37	12	

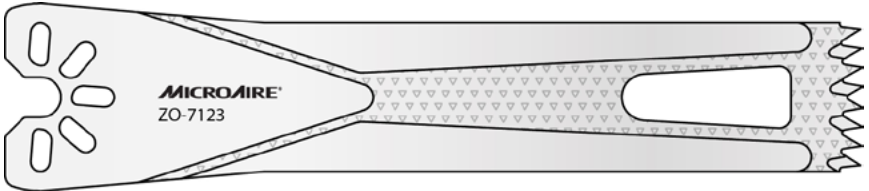
						MicroAire® and Hall® Large-Power Saws	StableCut®
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-72689	25	90	0.89	0.89	12		
ZO-72709	21	90	1.0	1.0	12		
ZO-72719	13	90	1.19	1.19	12		
ZO-72730	19	105	1.27	1.27	12		
ZO-72737	21	90	1.35	1.35	12		
ZO-73416	25	90	1.0	1.0	12		
ZO-73417	25	90	1.47	1.47	12		
ZO-73418	25	90	1.37	1.37	12		
ZO-73419	25	90	1.19	1.19	12		

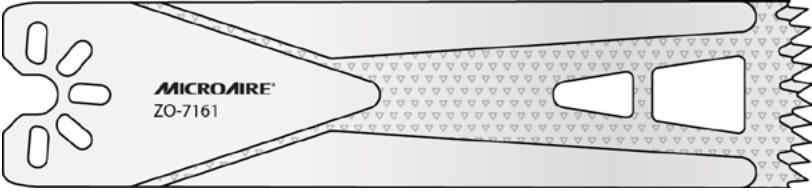
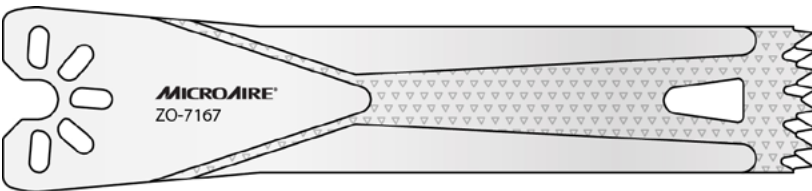
MicroAire® and Hall® Large-Power Saws StableCut®

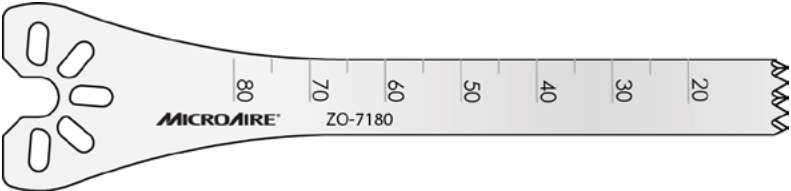
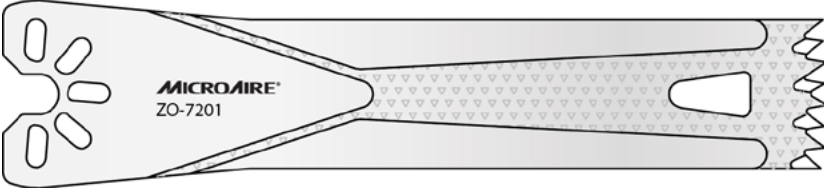
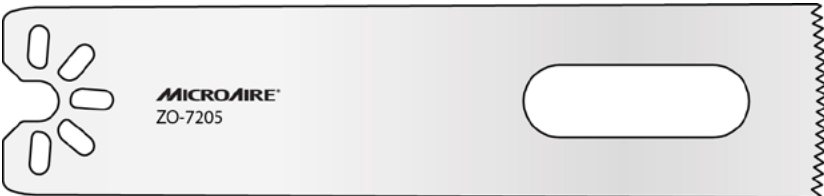
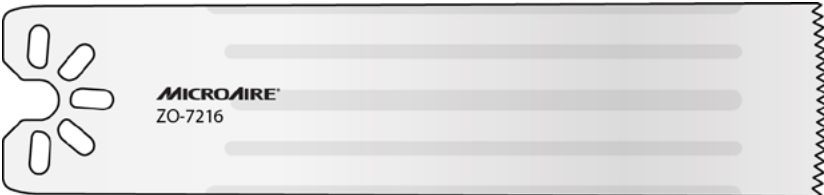
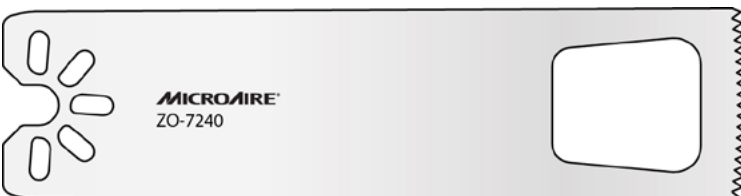
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-73420	13	90	1.27	1.27	12	
ZO-73421	13	90	1.0	1.0	12	
ZO-73422	13	90	1.37	1.37	12	
ZO-73489	13	90	0.89	0.89	12	
ZO-73613	21	90	1.19	1.19	12	
ZO-73616	21	90	1.27	1.27	12	

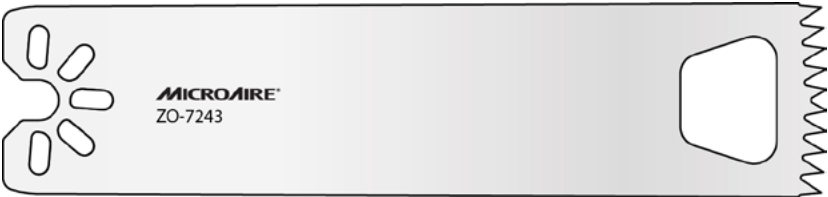
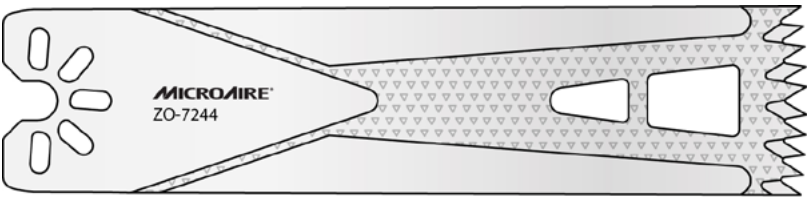
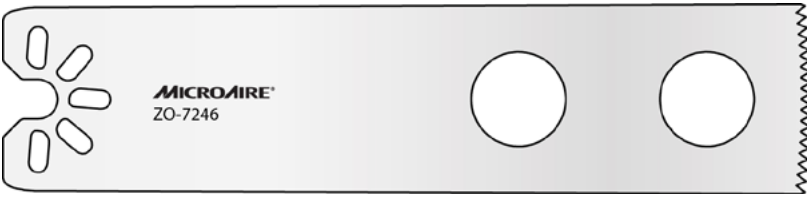
MicroAire® and Hall® Large-Power Saws All MicroAire® and Linvatec® (Except D-Cut®)

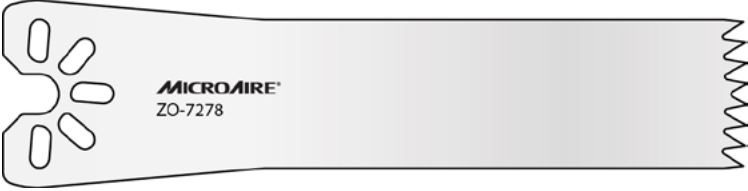
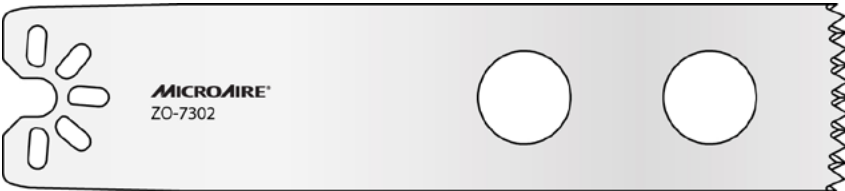
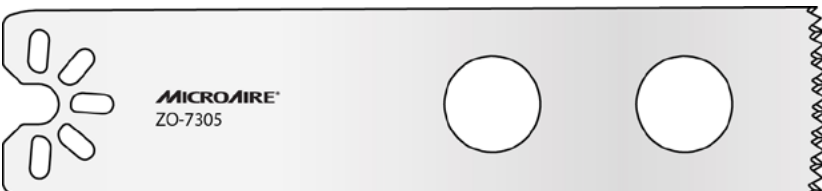
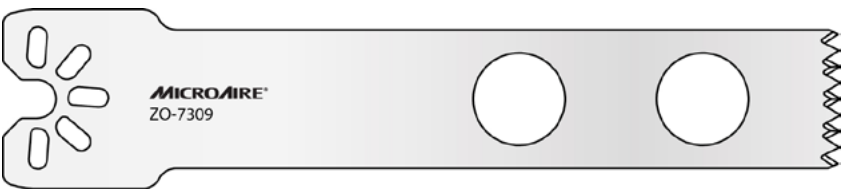
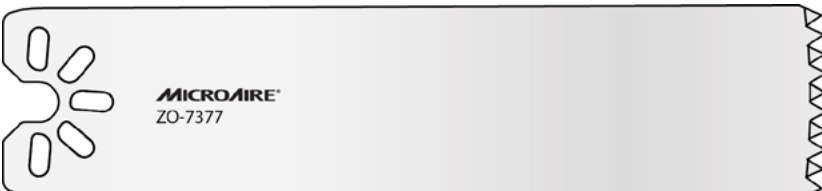
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-640	30.5	50	1.0	0.64	12	
ZO-7050	48	63	0.90	0.64	18	
ZO-7052	13	64	0.90	0.64	18	

MicroAire® and Hall® Large-Power Saws						All MicroAire® and Linvatec® (Except D-Cut®)
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-7053	19.5	63	0.90	0.64	18	
ZO-7055	32	68	1.50	1.20	18	
ZO-7111	25.4	90	1.40	1.27	12	
ZO-7123	19.5	95	1.27	1.27	10	
ZO-7145	11	30	1.10	0.64	18	
ZO-7151	33.3	63	0.90	0.64	18	

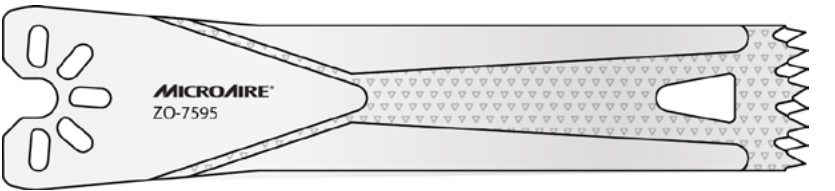
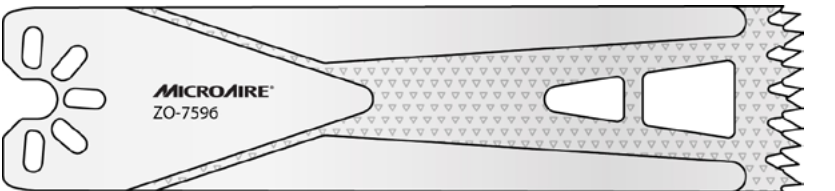
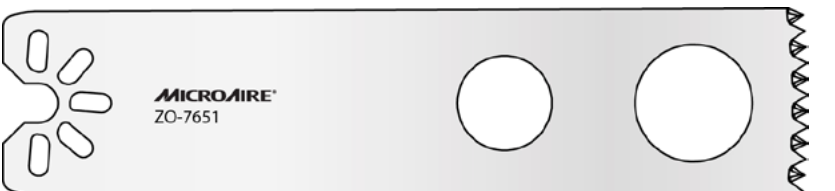
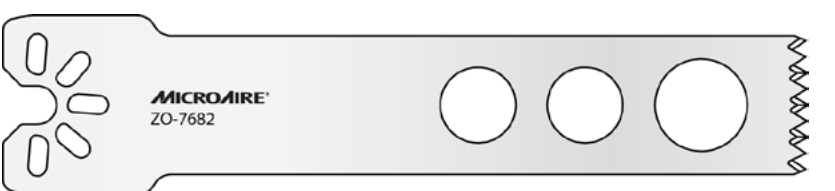
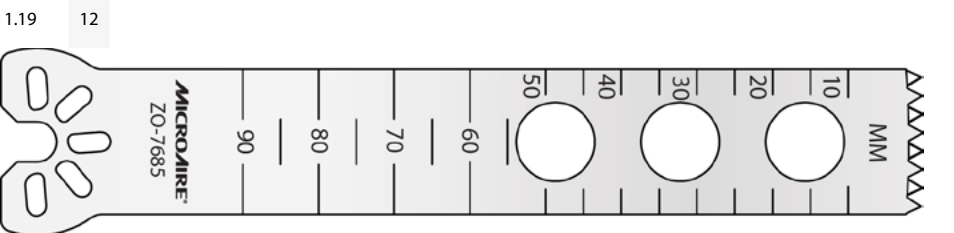
MicroAire® and Hall® Large-Power Saws						All MicroAire® and Linvatec® (Except D-Cut®)	
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-7154	19.5	71	1.35	1.20	18		
ZO-7156	33.5	64	0.90	0.80	18		
ZO-7161	25.4	90	1.27	1.27	10		
ZO-7163	25.4	90	1.37	1.37	10		
ZO-7165	25.4	90	1.47	1.47	10		
ZO-7166	32	63	0.90	0.80	10		
ZO-7167	19.5	90	1.20	1.20	10		
ZO-7172	19.5	71	0.90	0.80	18		

MicroAire® and Hall® Large-Power Saws						All MicroAire® and Linvatec® (Except D-Cut®)	
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-7180	10	85	1.27	1.27	12		
ZO-7201	19.5	90	1.27	1.27	10		
ZO-7205	25.4	90	1.0	0.90	12		
ZO-7216	25.4	90	1.35	1.35	12		
ZO-7240	25.4	80	1.40	1.20	12		
ZO-7241	19.5	80	1.40	1.27	10		
ZO-7242	25.4	90	1.40	1.27	10		

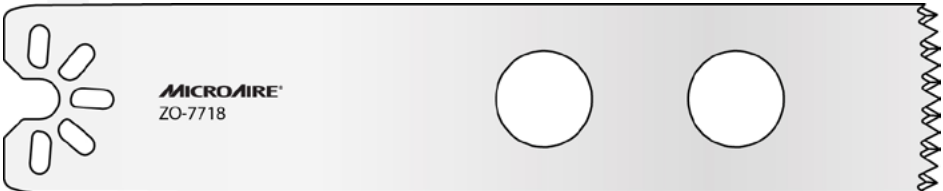
MicroAire® and Hall® Large-Power Saws						All MicroAire® and Linvatec® (Except D-Cut®)	
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-7243	25.4	90	0.90	0.80	10		
ZO-7244	25.4	90	1.50	1.50	10		
ZO-7245	13	70	1.50	1.27	12		
ZO-7246	25.4	90	0.90	0.90	12		
ZO-7250	13	90	1.0	0.90	16		
ZO-7251	12	89	1.35	1.35	12		
ZO-7274	19.5	86	0.90	0.80	18		

MicroAire® and Hall® Large-Power Saws						All MicroAire® and Linvatec® (Except D-Cut®)	
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-7277	19.5	71	0.90	0.80	10		
ZO-7278	19.5	80	0.9	0.80	10		
ZO-7302	25.4	95	0.89	0.89	12		
ZO-7305	25.4	90	1.20	1.20	12		
ZO-7309	19.1	90	0.89	0.89	12		
ZO-7377	25.4	90	1.50	1.50	12		
ZO-7380	19	85	1.35	1.35	12		

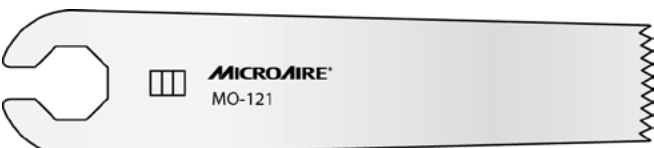
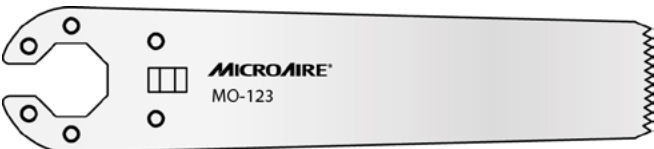
MicroAire® and Hall® Large-Power Saws All MicroAire® and Linvatec® (Except D-Cut®)

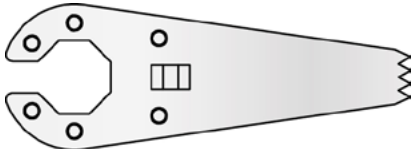
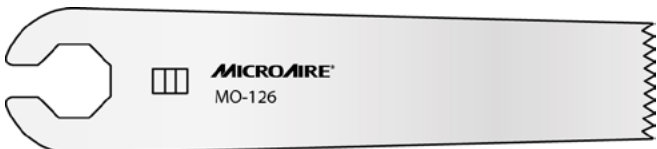
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-7390	19	90	1.35	1.35	12	
ZO-7595	19.5	90	0.90	0.90	10	
ZO-7596	25.4	90	0.90	0.90	10	
ZO-7651	25.4	90	1.19	1.19	12	
ZO-7654	13	90	1.19	1.19	12	
ZO-7682	19.1	85	1.27	1.27	12	
ZO-7685	19.1	105	1.19	1.19	12	

MicroAire® and Hall® Large-Power Saws All MicroAire® and Linvatec® (Except D-Cut®)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-7689	12.7	94	0.89	0.89	12	
ZO-7718	25	105	1.24	1.24	12	

3M® and Linvatec® Power Instruments (Models Maxi-Driver® G120A, L120A and Maxion™ E220)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
MO-121	15.2	50	1.60	1.0	10	
MO-122	31.5	50	1.60	1.0	10	
MO-123	15.2	50	1.10	0.60	13	

3M® and Linvatec® Power Instruments (Models Maxi-Driver® G120A, L120A and Maxion™ E220)						
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
MO-124	31.5	50	1.10	0.60	12	
MO-125	6.6	20	1.0	0.60	15	
MO-126	15.2	50	1.60	1.0	12	
MO-127	31.5	50	1.60	1.0	12	
MO-128	31.5	50	0.90	0.60	22	
MO-129	15.5	50	0.90	0.60	22	

Reciprocating Saw Blades

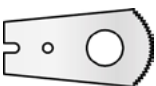
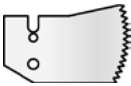
MicroAire® and Hall® Large-Power Saws

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZR-058	12.2	87.9	1.0	0.64	18	
ZR-059	9.4	52.8	1.0	0.64	18	
ZR-060	9.4	38.1	0.64	0.40	33	
ZR-061	12.2	80	1.80	1.27	18	
ZR-073	12.7	73.5	0.90	0.80	20	
ZR-079	12.7	76	0.90	0.90	16	
ZR-160	12.8	38.1	0.50	0.38	24	

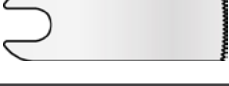
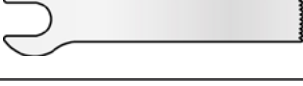
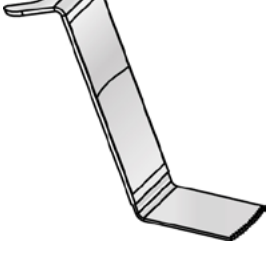
BLADES for Small-Bone Power Systems

Oscillating Saw Blades

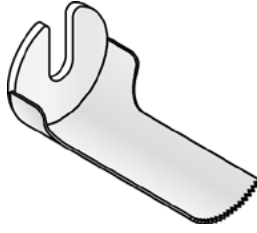
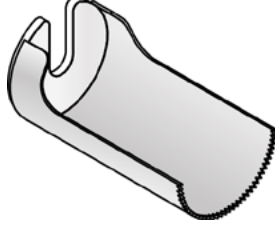
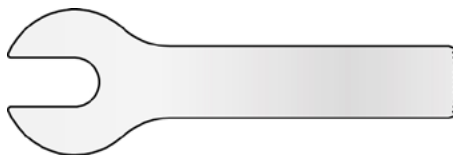
MicroAire® Oscillating Saw (Models 1700, 1710, 1970, 2700, 2710, 5922, 6700, and 6710)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	 CAUTION: Do not use the long 1200 Series sagittal blades in these models of oscillating saws, as they will overstress the oscillating mechanism.
1700-001	10.2	11	0.5	0.3	36	
1700-002	10.2	12	0.5	0.3	36	
1700-003	13	16.5	0.5	0.3	36	
1700-007	10.2	11	0.5	0.3	36	
1700-008	10.2	12	0.5	0.3	36	
1700-009	10.2	7.7	0.5	0.3	36	
1700-010	10.2	9.3	0.5	0.3	36	
1700-012	10.2	10.2	0.5	0.3	36	
1700-013	9.1	7	0.5	0.3	36	 

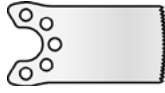
MicroAire® and Hall® Oscillating Saw (Models 1740, 1972, 2740, 5972 and 6740)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	 CAUTION: These blades cannot be used in MicroAire® oscillating saw models 1700, 1710, 1970, 1975, 1976, 2700, 2710, 6700, and 6710.
ZO-060	5.8	12	0.6	0.4	36	
ZO-061	9.7	12	0.6	0.4	36	
ZO-062	5.8	18.5	0.6	0.4	36	
ZO-063	9.7	18.5	0.6	0.4	36	
ZO-064	5.8	28.5	0.6	0.4	36	
ZO-065	9.7	28.5	0.6	0.4	36	
ZO-066	18.5	28.5	0.6	0.4	36	
ZO-067	5.8	13.5	0.6	0.4	36	
ZO-068	9.7	13.5	0.6	0.4	36	
ZO-069	9.7	13.5	0.6	0.6	36	

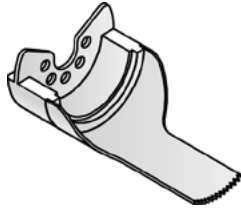
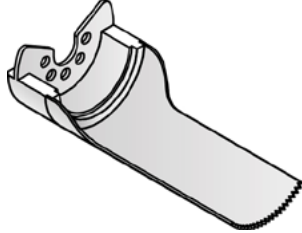
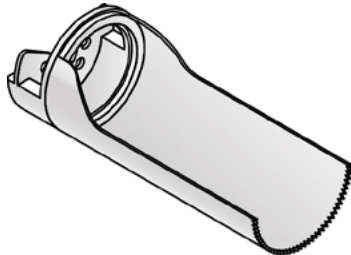
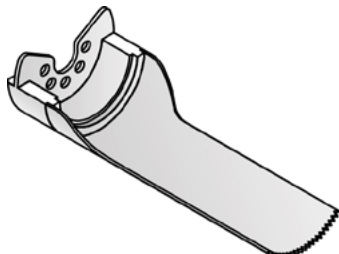
MicroAire® and Hall® Oscillating Saw (Models 1740, 1972, 2740, 5972 and 6740)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-070	9	11.9	0.6	0.4	33	
ZO-071	13.5	21.6	0.6	0.4	33	
ZO-072	18	22.9	0.6	0.4	33	
ZO-073	9.7	31.2	0.6	0.4	36	
ZO-076	45	16	0.7	0.4	18	
ZO-077	9.4	40.1	0.6	0.4	22	
ZO-165	9.7	27.6	0.6	0.4	36	

MicroAire® and Hall® Oscillating Saw (Models 1975 and 1976)

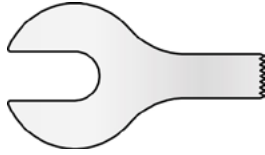
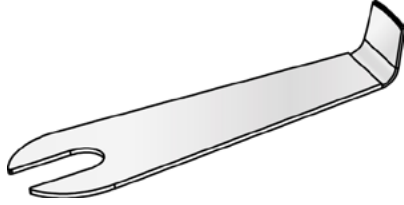
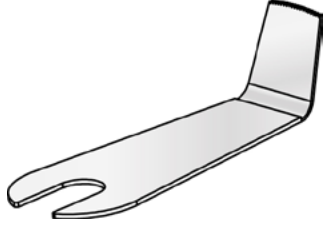
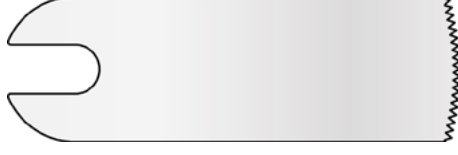
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-360	5.5	11	0.5	0.4	34	
ZS-361	10	11	0.5	0.4	34	
ZS-362	5.5	17.5	0.5	0.4	34	
ZS-363	10	17.5	0.5	0.4	34	
ZS-364	5.5	27.5	0.5	0.4	36	
ZS-365	10	27.5	0.5	0.4	36	
ZS-366	18.5	27	0.5	0.4	34	
ZS-367	5.8	13.5	0.6	0.4	36	
ZS-368	9.6	14	0.6	0.4	36	
ZS-369	10	33.5	0.5	0.4	34	

MicroAire® and Hall® Oscillating Saw (Models 1975 and 1976)

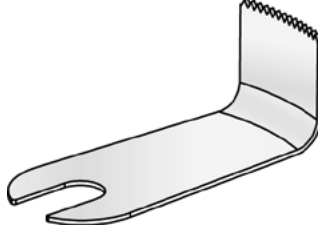
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-370	8.9	13.2	0.6	0.4	33	
ZS-371	13.5	23.1	0.6	0.4	33	
ZS-372	18	29	0.6	0.4	33	
ZS-373	13.5	28.7	0.6	0.4	33	

Sagittal Saw Blades

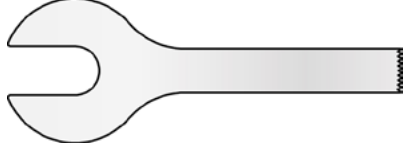
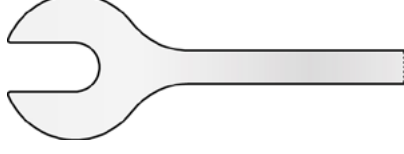
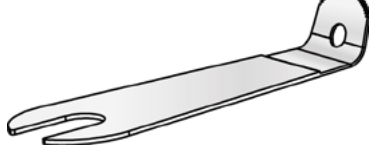
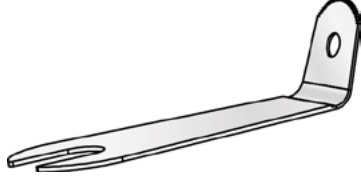
MicroAire® and Hall® Sagittal Saw (Models 1955, 2250, 5955, 6250, 6670, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-015	5.8	15	0.6	0.4	22	
ZS-018	9.7	24.4	0.6	0.4	36	
ZS-027	19.1	41	0.6	0.4	36	
ZS-028	9.4	7.1	0.6	0.4	36	
ZS-029	9.4	13.5	0.6	0.4	36	
ZS-030	19.1	25.4	0.6	0.4	22	
ZS-031	19.1	41	0.6	0.4	22	

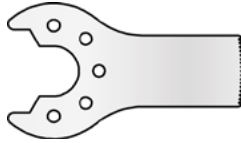
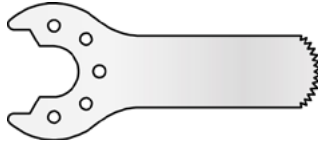
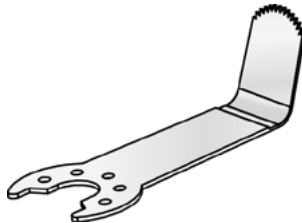
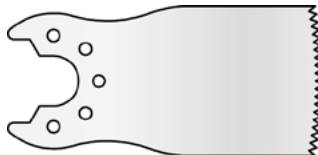
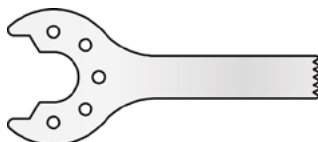
MicroAire® and Hall® Sagittal Saw (Models 1955, 2250, 5955, 6250, 6670, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-032	5.8	25.4	0.6	0.4	22	
ZS-033	14	41	0.6	0.4	22	
ZS-034	5.8	8.9	1.0	0.6	22	
ZS-035	5.8	17.3	1.0	0.6	22	
ZS-036	14	16	1.0	0.6	18	
ZS-037	14	25.4	0.6	0.4	22	
ZS-038	9.4	25.4	0.6	0.4	36	
ZS-038A	9.4	25.4	0.6	0.4	18	
XZS-038	9.5	25.4	1.0	0.8	22	

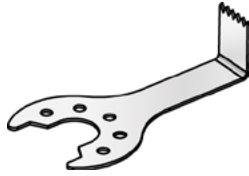
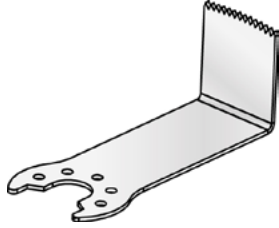
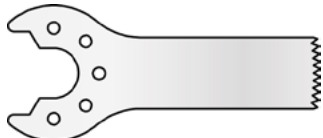
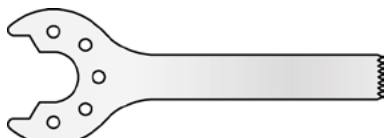
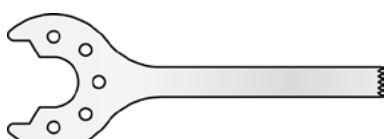
MicroAire® and Hall® Sagittal Saw (Models 1955, 2250, 5955, 6250, 6670, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-039	4.5	25.4	0.6	0.4	36	
ZS-040	9.4	41	0.6	0.4	36	
ZS-041	5.8	33.5	0.6	0.4	36	
ZS-042	4.3	35.6	0.6	0.4	36	
ZS-043	9.5	25.4	0.6	0.4	18	
ZS-104	9.4	15	0.6	0.4	36	
ZS-128	10.4	8.1	0.5	0.4	36	
ZS-129	10.4	13.5	0.5	0.4	36	

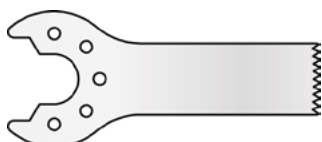
MicroAire® and Hall® Sagittal Saw (Models 1950, 1955, 2250, 5950, 5955, 6670, 6672, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-304	9.4	15	0.6	0.4	36	
ZS-318	9.5	25.4	0.5	0.4	31	
ZS-327	19.1	41	0.6	0.4	36	
ZS-329	9.7	14	0.56	0.4	34	
ZS-330	19.5	25.5	0.6	0.4	18	
ZS-331	19.5	41	0.6	0.4	18	
ZS-332	5.6	25.4	0.6	0.4	18	
ZS-333	14	41	0.6	0.4	18	

MicroAire® and Hall® Sagittal Saw (Models 1950, 1955, 2250, 5950, 5955, 6670, 6672, and 6673)

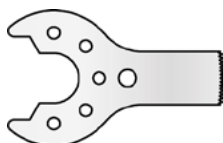
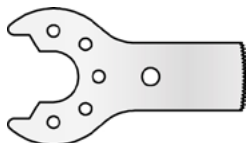
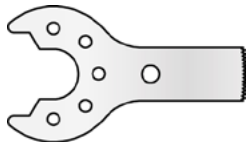
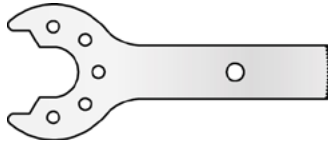
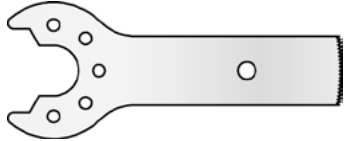
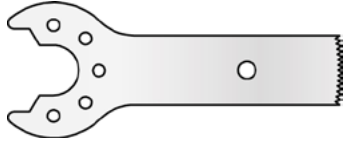
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-334	5.8	9	1.0	0.6	22	
ZS-336	14	16	1.0	0.6	18	
ZS-337	14	25.4	0.6	0.4	18	
ZS-338	9.4	25.4	0.5	0.4	34	
XZS-338	9.5	25.4	1.0	0.8	34	
ZS-339	4.5	25.4	0.6	0.4	34	
ZS-340	9.5	41	0.5	0.4	34	
ZS-341	6	34	0.5	0.4	34	
ZS-342	4	34	0.5	0.4	34	

MicroAire® and Hall® Sagittal Saw (Models 1950, 1955, 2250, 5950, 5955, 6670, 6672, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-343	9.5	25.4	0.6	0.4	18	

MicroAire® TPS, CORE, RemB, and Command 2 Stryker-style Blades

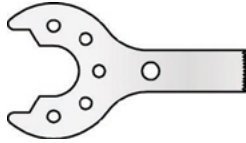
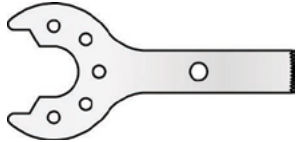
MicroAire® and Hall® Sagittal Saw (Models 1950, 1955, 2250, 5950, 5955, 6670, 6672, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-103	7	13	0.4	0.4	32	
ZS-105	9	17	0.4	0.4	32	
ZS-114	7	16	0.4	0.4	32	
ZS-115	7	29.5	0.4	0.4	32	
ZS-125	9	31	0.4	0.4	32	
ZS-225	9	31	0.4	0.4	22	
ZS-410	5.5	9.5	0.4	0.4	32	

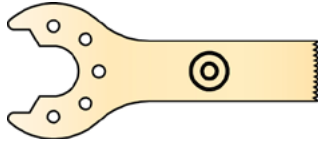
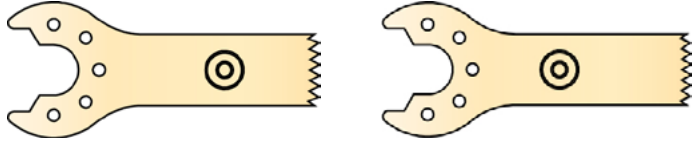
MicroAire® TPS, CORE, RemB, and
Command 2 Stryker-style Blades

MicroAire® and Hall® Sagittal Saw

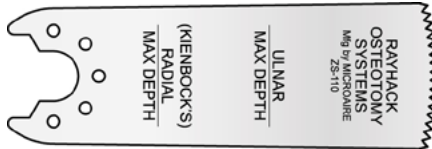
(Models 1950, 1955, 2250, 5950, 5955,
6670, 6672, and 6673)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-412	5.5	18	0.4	0.4	32	
ZS-414	5.5	22	0.4	0.4	32	

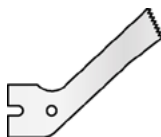
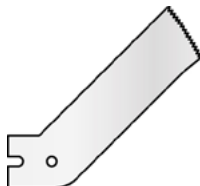
Titanium Nitride-Coated ACL Blades

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-382	8	12	0.6	0.4	18	
ZS-390	9.5	10	0.6	0.4	18	
ZS-395	9.5	15	0.6	0.4	18	

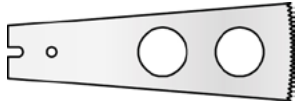
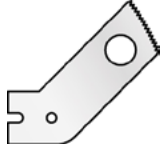
For the Wright® Medical Technologies Rayhack® System

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZS-110	19.1	40	0.6	0.5	22	
ZS-510	14	40	0.6	0.5	22	

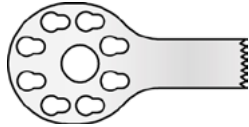
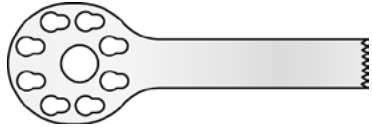
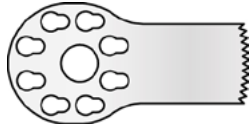
MicroAire® Sagittal Saw (Models 1220, 1920T, 1920, 1922, 2200, 5922 and 6200)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
1200-001	12.7	19.8	0.5	0.3	36	
1200-003	12.7	24.3	0.5	0.3	36	
1200-005	6.9	19.8	0.5	0.3	36	
1200-006	4.1	17	0.5	0.3	36	
1200-007	4.3	15	0.5	0.3	32	
1200-008	7.9	21.5	0.5	0.3	36	
1200-009	10.2	17.8	0.5	0.3	36	
1200-010	10.2	13.7	0.5	0.3	36	
1200-011	8.1	6.6	0.5	0.3	36	
1200-012	9.1	13.7	0.5	0.3	36	

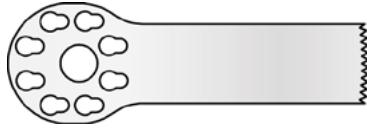
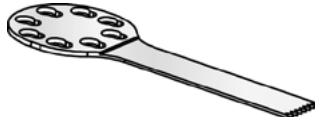
MicroAire® Sagittal Saw (Models 1220, 1920T, 1920, 1922, 2200, 5922 and 6200)

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
1200-013	12.7	29	0.5	0.3	36	
1200-014	7.5	15	0.5	0.3	36	
1200-017	5.8	18	0.5	0.3	36	
1200-018	10.2	18.3	0.5	0.3	36	
1200-019	6.9	28.7	0.5	0.3	36	
1200-019A	6.9	13.5	0.5	0.3	22	
1200-038	9.4	25.5	0.5	0.3	36	

Hall®(3M®) Mini-Driver K120 Sagittal Saw Insert

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
MS-121	6.4	15.7	1.0	0.7	27	
MS-122	6.4	31.7	1.0	0.7	26	
MS-123	10.4	15.7	1.0	0.7	27	

Hall®(3M®) Mini-Driver K120 Sagittal Saw Insert

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
MS-124	10.4	31.7	1.0	0.7	26	
MS-125	15.5	15.7	1.0	0.7	27	
MS-131	6.4	15.7	0.7	0.7	27	
MS-132	6.4	31.7	0.7	0.7	26	
MS-134	10.4	31.7	0.7	0.7	26	

Reciprocating Saw Blades

MicroAire® Reciprocating Saw Models 1400, 1945, 2420, 5945, 6400 and 6420

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
1400-001	8	25.4	0.7	0.4	32	
1400-002	8	38.1	0.7	0.4	32	
1400-003	4.3	12.7	0.9	0.5	30	
1400-004	4.3	12.7	0.9	0.5	30	
1400-005	5	22.9	0.6	0.4	30	
1400-006	9	40	0.8	0.5	24	

MicroAire® Reciprocating Saw Models 1945, 2420, 5945, 6400 and 6420						
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZR-020	7.6	25.4	0.6	0.4	33	
ZR-021	9.4	39.1	0.6	0.4	33	
ZR-050	5.8	6.8	0.8	0.5	21	
ZR-051	5.0	17.2	0.8	0.5	21	
ZR-052	6.1	23.1	0.8	0.5	21	
ZR-053	8.1	25.1	0.8	0.5	21	
ZR-120	8.1	25.4	0.5	0.4	22	
ZR-154	7.4	23.6	0.6	0.4	36	
ZR-733	6.4	22.9	0.6	0.4	20	

STERNUM SAW BLADES

Sternum Saw Blades

Linvatec®

Sternum Saw

Model 5059-05

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZR-031	9.3	53.6	1.0	0.64	18	

MicroAire®

Sternum Saw

Model 4800

Linvatec®

Sternum Saw

Model 5059-05

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZR-032	9.4	34	1.0	0.6	18	

MicroAire®

Sternum Saw

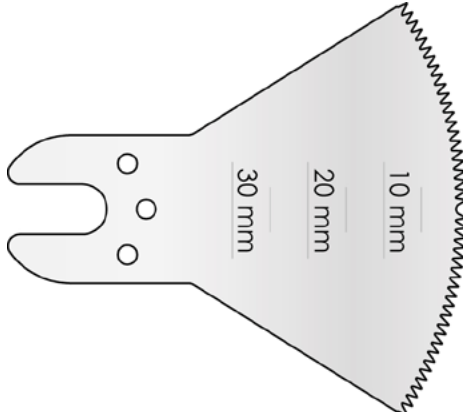
Model 7800

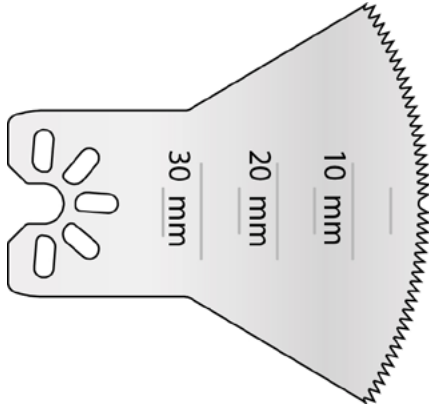
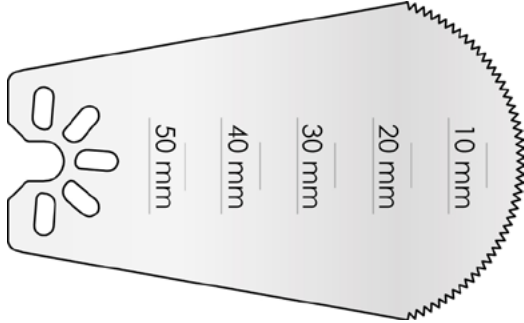
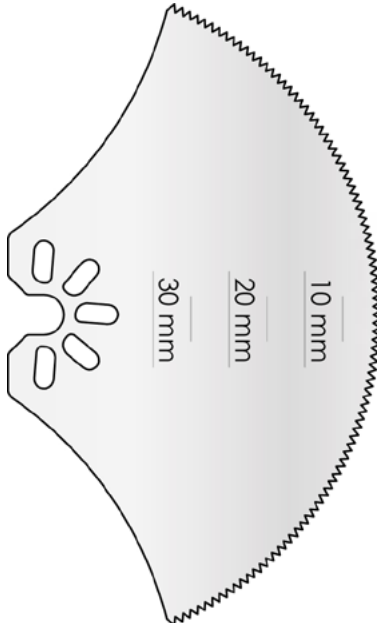
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZR-032M	9.4	34	1.0	0.6	18	
ZR-033M	9.4	34	1.0	0.6	16	

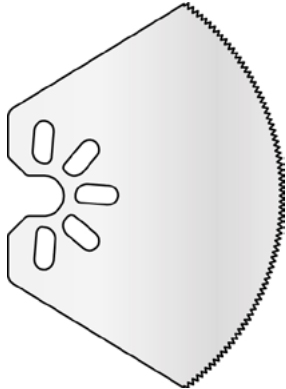
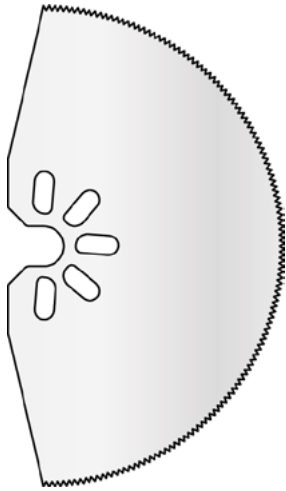
Linvatec® / Hall®

Repeat Sternotomy Saw

Model 5057

Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI	
ZO-079	53.3	30	0.6	0.4	18	

MicroAire® Oscillating Saw		Models 7200, 7205, 7600 and 7506		Linvatec® Oscillating Saw		Models 5052 and 5044-02 PowerPro®, M-Power®, and M-Power®2	
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-7079	53	35	0.6	0.4	22		
ZO-7125	42	50	1.1	0.64	22		
ZO-7137	83	30	1.0	0.64	22		

MicroAire® Oscillating Saw		Models 7200, 7205, 7600 and 7506		Linvatec® Oscillating Saw		Models 5052 and 5044-02 PowerPro®, M-Power®, and M-Power®2	
Part No.	Width (mm)	Usable Length (mm)	Cut Thickness (mm)	Material Thickness (mm)	TPI		
ZO-7225	51	19	0.7	0.4	18		
ZO-7226	63.5	19	0.7	0.4	18		

Medium-Length Burs

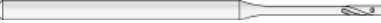
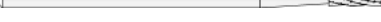
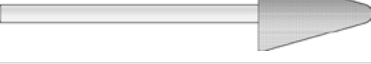
MicroAire®

Drill with Built-on
Medium Bur Guard

Models 1140, 1910, 1910T, 1930, 1932, 1940, 1940T, 2100, 2120, 2130, 2910, 2950, 5930, 6100, 6120, and 6130. Also fits Linvatec® Micro Drills with Medium Bur Guard.

Type	Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	 CAUTION: Do not use burs with head diameters larger than 4mm in the 1910, 1910T, 2130, 2910, 2950, or 6130 High Speed Drills. The burs may break and cause injury. Length must be between 47.8mm–56.9mm. Do not use long burs without the proper bur guard.
End-Side Cut Carbide (Suture Holes)	ZB-100	0.8	4.6	6	50.8	
Side-Cut Carbide	ZB-101	1	5.1	6	50.8	
	ZB-104	1.5	7.6	6	50.8	
Side-Cut Tapered Carbide	ZB-102	1.5	10.2	6	50.8	
	ZB-103	1	5.1	6	50.8	
	ZB-106	2	7.6	6	50.8	
	ZB-113	2	5	6	50.8	
	ZB-114	2	19.1	6	56.9	
Lindemann Stainless Steel	ZB-112	2.3	31.5	4	66.5	
Round-Cut Carbide	ZB-121	1	1	6	50.8	
	ZB-122	1.5	1.5	6	50.8	
	ZB-124	2	2	6	50.8	
	ZB-125	2.4	2.4	6	50.8	
	ZB-126	3	3	8	50.8	
	ZB-128	4	4	8	48.3	
	ZB-130	5	5	10	48.3	
	ZB-132	6.5	6.5	10	48.3	
Oval-Cut Carbide	ZB-136	4	7.9	8	47.8	
	ZB-138	5.5	9.9	10	47.8	
Swanson Pilot Point Reamer	ZB-142	2	8.6	4	57.2	
	ZB-143	3	13.7	5	61.2	
	ZB-144	4	16.8	6	62	
	ZB-145	5	18.5	8	63	

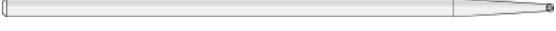
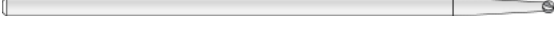
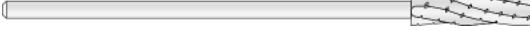
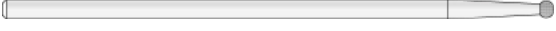
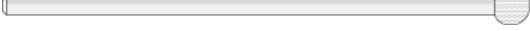
MicroAire® Drill with Built-on Medium Bur Guard Models 1140, 1910, 1910T, 1930, 1932, 1940, 1940T, 2100, 2120, 2130, 2910, 2950, 5930, 6100, 6120, and 6130. Also fits Linvatec® Micro Drills with Medium Bur Guard.

Type	Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	 CAUTION: Do not use burs with head diameters larger than 4mm in the 1910, 1910T, 2130, 2910, 2950, or 6130 High Speed Drills. The burs may break and cause injury. Length must be between 47.8mm–56.9mm. Do not use long burs without the proper bur guard.
Wire Pass Drill	ZB-147	1.1	19.1	2	50.8	
	ZB-148	1.5	19.1	2	50.8	
	ZB-149	2	19.1	2	50.8	
	ZB-447	1.1	19.1	2	36.8	
	ZB-448	1.5	19.1	2	36.8	
	ZB-449	2.0	19.1	2	36.8	
Barrel-Cut Carbide	ZB-152	2	10.2	8	50.8	
	ZB-154	4	10.2	12	47.8	
Round Diamond	ZB-160	2	2	N/A	50.8	
	ZB-162	3	3	N/A	50.8	
	ZB-164	4	4	N/A	47.8	
	ZB-168	6	6	N/A	47.8	
Oval Diamond	ZB-170	4	7.9	N/A	48.3	
	ZB-180	4	15	N/A	48.3	
	ZB-181	7	15	N/A	48.3	
Diamond Dermabraders	ZB-182	10.5	15	N/A	49	
	ZB-183	10.2	15	N/A	49.3	
	ZB-184	5	10.2	N/A	49.3	
Lindemann Carbide	ZB-190	1.9	15.2	4	48.3	
Shannon 44 (Stainless Steel)	ZB-192	2.3	14	3	50.8	

Long Burs

Drill with 1100-005 Bur Guard 1140, 1930, 1932, 1940, 1940T, 2100, 2120, 5930 6100, and 6120
 Drill with 2910-011 Bur Guard 2130, 2910, 2950, and 6130 high-speed drills
 Drill with 1910-011 Bur Guard 1910 and 1910T high-speed drills
 Also fits Linvatec® Micro Drills with Long Bur Guard.

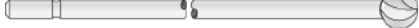
MicroAire®

All require long bur guard.						 CAUTION: Do not use burs with head diameters larger than 4mm in the 1910, 1910T, 2130, 2910, 2950, or 6130 High Speed Drills. The burs may break and cause injury. Length must be between 69.8mm–80mm. Do not use long burs without the proper bur guard.
Type	Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
Lindemann Carbide	ZB-200	1.9	15.2	4	69.9	
Side-Cut Tapered Carbide	ZB-203	1.7	12.2	6	69.9	
Side-Cut Carbide	ZB-204	1.5	7.6	6	73.2	
Side-Cut Tapered Carbide	ZB-206	2	7.6	6	73.2	
	ZB-214	2	19.1	6	80	
Round-Cut Carbide	ZB-220	1	1	6	73.2	
	ZB-222	1.5	1.5	6	73.2	
	ZB-224	2	2	6	73.2	
	ZB-226	3	3	8	73.2	
	ZB-228	4	4	8	69.9	
	ZB-230	5	5	10	69.9	
	ZB-232	6.5	6.5	10	69.9	
Oval-Cut Carbide	ZB-236	4	7.9	8	69.9	
	ZB-238	5.5	9.9	10	69.9	
Barrel-Cut Carbide	ZB-252	2	10.2	8	73.2	
	ZB-254	4	10.2	12	69.9	
Round Diamond	ZB-260	2	2	N/A	73.2	
	ZB-262	3	3	N/A	73.2	
	ZB-264	4	4	N/A	69.1	
	ZB-266	5	5	N/A	69.1	
	ZB-268	6	6	N/A	69.1	
Oval Diamond	ZB-270	4	7.9	N/A	69.9	
Shannon 44 (Stainless Steel)	ZB-292	2.3	14	3	69.9	

Extra-Long Burs

Drill with 1100-005 Bur Guard 1140, 1930, 1932, 1940, 1940T, 2100, 2120, 5930 6100, and 6120
 Drill with 2910-019 Bur Guard 2130, 2910, 2950, and 6130 high-speed drills
 Drill with 1910-019 Bur Guard 1910 and 1910T high-speed drills
 Also fits Linvatec® Micro Drills with Extra-Long Bur Guard.

MicroAire®

Type	Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
CAUTION: Do not use burs with head diameters larger than 4mm in the 1910, 1910T, 2130, 2910, 2950, or 6130 High Speed Drills. The burs may break and cause injury. Do not use long burs without the proper bur guard.						
Side-Cut Tapered Carbide	ZB-306	2	7.6	6	95.3	
Round-Cut Carbide	ZB-322	1.5	1.5	6	95.3	
	ZB-324	2	2	6	95.3	
	ZB-326	3	3	8	95.3	
	ZB-328	4	4	8	95.3	
	ZB-330	5	5	10	95.3	
Oval-Cut Carbide	ZB-336	4	7.9	8	95.3	
	ZB-338	5.5	9.9	10	95.3	
Round Diamond	ZB-360	2	2	N/A	94.5	
	ZB-362	3	3	N/A	94.5	
	ZB-364	4	4	N/A	95.3	
	ZB-366	5	5	N/A	94.5	

Oral Surgery Burs

MicroAire® Drill

Models 1140, 1910, 1910T, 1930, 1932, 1940, 1940T, 2100, 2120, 2130, 2910, 2950, 5930, 6100, 6120, and 6130 high-speed drills. Also fits Linvatec® Micro and High Speed Drills with Medium Bur Guard.

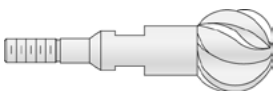
Type	Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
Round-Cut Carbide	ZB-70	2.3	2.3	8	50.8	
Side-Cut Carbide	ZB-71	1.6	4.4	6	50.8	
Tapered Side-Cut Carbide	ZB-72	1.2	4	6	50.8	
	ZB-73	1.6	4.6	6	50.8	
	ZB-74	2.1	4.6	6	50.8	
Round-Cut Carbide	ZB-75	1.8	1.6	6	50.8	
Round-End Carbide	ZB-76	1.2	4.3	6	50.8	
Side-Cut Carbide	ZB-77	1.2	4.3	6	50.8	
Cross-Cut Fissure	ZB-78	2.1	4.6	8	50.8	

*See page A-39 for Drill Nose Shield

Cebotome™ Drill Burs

Linvatec® (Hall®)

Cebotome™ drills

Type	Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
Side Cutting	ZB-582	6	24.3	2	42.7	
Round Cutting	ZB-584	7	6.1	10	30.2	
	ZB-585	8	8	10	31.3	
	ZB-586	9	8.4	10	32.5	

Short Round Burs

Linvatec® (Hall®)

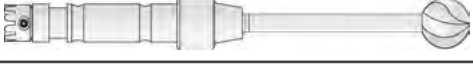
Osteon™ drill system

Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
ZB-617	6.5	6.5	10	30	
ZB-618	5	5	10	30	
ZB-619	4	4	8	30	

Medium Round Burs

Linvatec® (Hall®)

Osteon™ drill system

Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
ZB-625	6.5	6.5	10	38	
ZB-626	5	5	10	38	
ZB-627	4	4	8	38	
ZB-628	3	3	8	38	
ZB-647	5	5	0	38	
ZB-648	4	4	0	38	

Long Round Burs

Linvatec® (Hall®)

Osteon™ drill system

Part No.	Head Size Diameter (mm)	Head Length (mm)	# of Flutes	Overall Length (mm)	
ZB-640	1.5	1.5	6	48	
ZB-642	2	2	6	48	
ZB-644	5	5	10	48	
ZB-645	4	4	8	48	
ZB-646	3	3	8	48	
ZB-650	5	5	0	48	
ZB-659	3	3	0	48	
ZB-662	2	2	0	48	
ZB-664	1	1	6	48	



1100-008

MicroAire® Drill Nose Shield

Special autoclavable silicone nose cap prevents blood and saline from entering the drill nose. Especially important during oral surgery where deposits can damage the bearing and gripping mechanism, resulting in overheating.

- For models 1930, 1932, 1940, 1940T, 2100, 2120, 2130, 2910, 2950, 5130, 5930, 5932, 6100, 6120, and 6130 drills
- Package of 10

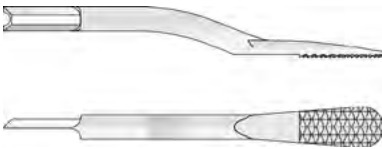
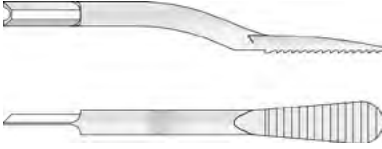
RASPS

Rasps

MicroAire®

Reciprocating Saw

Model 1400

Type	Part No.	Width (mm)	Usable Length (mm)	
Cross-Serrated	1400-011	5.3	11	
Cottle	1400-013	5.6	15.5	

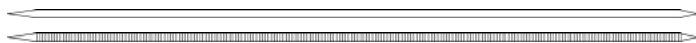
MicroAire®

Reciprocating Saw

Models 1945, 2420, 1915, 1915T, 5945, 6400 and 6420

Type	Part No.	Width (mm)	Usable Length (mm)	
Small Straight-Cut	ZR-080	5.5	9.5	
Large Straight-Cut	ZR-081	6.5	12.5	
Small Cross-Cut	ZR-082	5.5	9.5	
Large Cross-Cut	ZR-083	6.5	12.5	
Small Intramedullary	ZR-084	5.3	24.5	
Large Intramedullary	ZR-085	9.5	30.5	

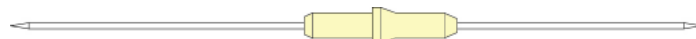
K-WIRES + STEINMANN PINS



K-Wires

6-packs only. To order non-sterile, add "NS" to the end of the product number (e.g. 1600-XXXNS)

Length	Point Style	Type	Wire Diameter				
			0.7 mm / .028"	0.9 mm / .035"	1.1 mm / .045"	1.4 mm / .054"	1.6 mm / .062"
4"	Single Trocar	Smooth	1604-128	1604-135	1604-145	1604-154	1604-162
	Double Trocar	Smooth	1600-428	1600-435	1600-445	1600-454	1600-462
		Threaded		1600-435T	1600-445T	1600-454T	1600-462T
6"	Double Trocar	Smooth	1600-628	1600-635	1600-645	1600-654	1600-662
		Threaded		1600-635T	1600-645T	1600-654T	1600-662T
9"	Single Trocar	Smooth		1600-9355	1600-9455	1600-9545	1600-9625
	Double Trocar	Smooth		1600-935	1600-945	1600-954	1600-962
		Threaded		1600-935T	1600-945T	1600-954T	1600-962T



K-Wires with Wire Guide

MicroAire® Pencil Grip Wire Driver Models 1620, 1625, 1995 and 2600
6-packs only

Length	Point Style	Type	Wire Diameter				
			0.7 mm / .028"	0.9 mm / .035"	1.1 mm / .045"	1.4 mm / .054"	1.6 mm / .062"
5.5"	Double Trocar	Smooth	1600-020	1600-021	1600-022		1600-023
			1600-001NS	1600-002NS	1600-003NS		1600-004NS
 Wire Guide only (non-sterile) Color			1600-009NS Red	1600-010NS Blue	1600-011NS Yellow	1600-013NS White	1600-012NS Black



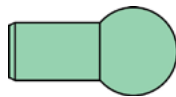
Steinmann Pins

6-packs only. To order non-sterile, add "NS" to the end of the product number (e.g. 1600-XXXNS)

Length	Point Style	Type	Wire Diameter						
			2.0 mm / 5/64"	2.4 mm / 3/32"	2.8 mm / 7/64"	3.2 mm / 1/8"	3.6 mm / 9/64"	4.0 mm / 5/32"	4.8 mm / 3/16"
9"	Single Trocar	Smooth	1620-509	1624-509	1628-509	1632-509	1636-509	1640-509	1648-509
		Threaded	1620-509T	1624-509T	1628-509T	1632-509T	1636-509T	1640-509T	1648-509T
	Double Trocar	Smooth	1620-109	1624-109	1628-109	1632-109	1636-109	1640-109	1648-109
		Threaded	1620-109T	1624-109T	1628-109T	1632-109T	1636-109T	1640-109T	1648-109T
	Single Diamond	Smooth	1620-609	1624-609	1628-609	1632-609	1636-609	1640-609	1648-609
		Threaded	1620-609T	1624-609T	1628-609T	1632-609T	1636-609T	1640-609T	1648-609T
	Double Diamond	Smooth	1620-709	1624-709	1628-709	1632-709	1636-709	1640-709	1648-709
		Threaded	1620-709T	1624-709T	1628-709T	1632-709T	1636-709T	1640-709T	1648-709T

K-Wire + Pin Protectors

Molded thermoplastic, protection from sharp points, and sold in sterile 4-packs.

For K-Wires	Wire Diameter						
	0.7 mm / .028"	0.9 mm / .035"	1.1 mm / .045"	1.6 mm / .062"			
	1601-028	1601-035	1601-045	1601-062			
Color	Yellow	Blue	White	Green			
For Steinmann Pins	Wire Diameter						
	2.0 mm / 5/64"	2.4 mm / 3/32"	2.8 mm / 7/64"	3.2 mm / 1/8"	3.6 mm / 9/64"	4.0 mm / 5/32"	4.8 mm / 3/16"
	1601-564	1601-332	1601-764	1601-180	1601-964	1601-532	1601-316
Color	Pink	Gray	Tan	Lavendar	Brown	Navy Blue	Orange

TWIST DRILLS + TAPS

Standard Orthopaedic Twist Drills

MicroAire® Jacobs® Drill Couplers		Also fits all other standard Jacobs®-style Drill Couplers (Hall®, Stryker®, Synthes®, Aesculap®, DeSouter®, Brasseler®, etc.)			To order non-sterile, add "NS" to the end of the product number (e.g. 8054-XXXNS)	
Type	Diameter (mm)	Length (mm)	MicroAire Part #	Replaces Zimmer®		
MicroAire Standard Twist Drills	1.6	127	8054-010	256-01-21		
	2.0		8054-011	256-01-22		
	2.4		8054-012	256-01-23		
	2.8		8054-013	256-01-24		
	3.2		8054-014	256-01-25		
	3.6		8054-015	256-01-26		
	4.0	127	8054-016	256-01-27		
	4.5		8054-017	256-01-28		
	4.8		8054-018	256-01-29		
	6.4		8054-019			
ECT® Twist Drills*	2.0	127	8054-020	2318-20		
	2.8		8054-021	2318-27		
	3.2		8054-022	2318-32		
	3.6		8054-023	2318-35		
	4.5		8054-024	2318-45		

*ECT® European Compression Technique

Small Power Instrument Twist Drills

MicroAire® Wire Driver			Models 1620, 1625, 1925, 1995 and 2600 with wire guide	Wire guide required (not included). Sold individually. Sterile or non-sterile; to order non-sterile, add "NS" to the end of the product number (e.g. 8051-XXXNS).
Diameter (mm)	Length (mm)	MicroAire Part #		
2.0	65	8051-120		
2.7	79	8051-127		
Wire Guide for Drills Sold in packs of 6 Non-sterile only		1600-012NS		



Standard Twist Drills with 2.34 mm Dental Bur Shank

May be used in standard MicroAire® and Linvatec® Micro Drills

Sold individually. Sterile or non-sterile; to order non-sterile, add "NS" to the end of the product number (e.g. 8051-XXXNS)

Diameter (mm)	Length (mm)	MicroAire Part #	CAUTION: Do not operate above 10,000 RPM. Do not use in high-speed drill.
1.1	49	8051-011	
1.3	49	8051-013	
1.5	49	8051-015	
1.8	49	8051-018	
2.0	65	8051-020	
2.2	49	8051-021	
2.7	79	8051-027	
3.2	79	8051-032	
3.5	79	8051-035	

INDEX BY PART NUMBER

Part #	Page	Part #	Page	Part #	Page	Part #	Page	Part #	Page
1200-001	A-27	1600-002NS	A-41	1600-645T	A-41	1601-764	A-42	1628-509T	A-42
1200-003	A-27	1600-003NS	A-41	1600-654	A-41	1601-964	A-42	1628-609	A-42
1200-005	A-27	1600-004NS	A-41	1600-654T	A-41	1604-128	A-41	1628-609T	A-42
1200-006	A-27	1600-009NS	A-41	1600-662	A-41	1604-135	A-41	1628-709	A-42
1200-007	A-27	1600-010NS	A-41	1600-662T	A-41	1604-145	A-41	1628-709T	A-42
1200-008	A-27	1600-011NS	A-41	1600-935	A-41	1604-154	A-41	1632-109	A-42
1200-009	A-27	1600-012NS	A-41	1600-935T	A-41	1604-162	A-41	1632-109T	A-42
1200-010	A-27	1600-012NS	A-44	1600-945	A-41	1620-109	A-42	1632-509	A-42
1200-011	A-27	1600-013NS	A-41	1600-945T	A-41	1620-109T	A-42	1632-509T	A-42
1200-012	A-27	1600-020	A-41	1600-954	A-41	1620-509	A-42	1632-609	A-42
1200-013	A-28	1600-021	A-41	1600-954T	A-41	1620-509T	A-42	1632-609T	A-42
1200-014	A-28	1600-022	A-41	1600-962	A-41	1620-609	A-42	1632-709	A-42
1200-017	A-28	1600-023	A-41	1600-962T	A-41	1620-609T	A-42	1632-709T	A-42
1200-018	A-28	1600-428	A-41	1600-9355	A-41	1620-709	A-42	1636-109	A-42
1200-019	A-28	1600-435	A-41	1600-9455	A-41	1620-709T	A-42	1636-109T	A-42
1200-019A	A-28	1600-435T	A-41	1600-9545	A-41	1624-109	A-42	1636-509	A-42
1200-038	A-28	1600-445	A-41	1600-9625	A-41	1624-109T	A-42	1636-509T	A-42
1400-001	A-29	1600-445T	A-41	1601-028	A-42	1624-509	A-42	1636-609	A-42
1400-002	A-29	1600-454	A-41	1601-035	A-42	1624-509T	A-42	1636-609T	A-42
1400-003	A-29	1600-454T	A-41	1601-045	A-42	1624-609	A-42	1636-709	A-42
1400-004	A-29	1600-462	A-41	1601-062	A-42	1624-609T	A-42	1636-709T	A-42
1400-005	A-29	1600-462T	A-41	1601-180	A-42	1624-709	A-42	1640-109	A-42
1400-006	A-29	1600-628	A-41	1601-316	A-42	1624-709T	A-42	1640-109T	A-42
1400-011	A-40	1600-635	A-41	1601-332	A-42	1628-109	A-42	1640-509	A-42
1400-013	A-40	1600-635T	A-41	1601-532	A-42	1628-109T	A-42	1640-509T	A-42
1600-001NS	A-41	1600-645	A-41	1601-564	A-42	1628-509	A-42	1640-609	A-42

Part #	Page	Part #	Page	Part #	Page	Part #	Page	Part #	Page
1640-609T	A-42	8051-035	A-44	MS-122	A-28	ZB-122	A-34	ZB-190	A-35
1640-709	A-42	8051-120	A-44	MS-123	A-28	ZB-124	A-34	ZB-192	A-35
1640-709T	A-42	8051-127	A-44	MS-124	A-29	ZB-125	A-34	ZB-200	A-36
1648-109	A-42	8054-010	A-43	MS-125	A-29	ZB-126	A-34	ZB-203	A-36
1648-109T	A-42	8054-011	A-43	MS-131	A-29	ZB-128	A-34	ZB-204	A-36
1648-509	A-42	8054-012	A-43	MS-132	A-29	ZB-130	A-34	ZB-206	A-36
1648-509T	A-42	8054-013	A-43	MS-134	A-29	ZB-132	A-34	ZB-214	A-36
1648-609	A-42	8054-014	A-43	XZS-038	A-21	ZB-136	A-34	ZB-220	A-36
1648-609T	A-42	8054-015	A-43	XZS-338	A-24	ZB-138	A-34	ZB-222	A-36
1648-709	A-42	8054-016	A-43	ZB-70	A-37	ZB-142	A-35	ZB-224	A-36
1648-709T	A-42	8054-017	A-43	ZB-71	A-37	ZB-143	A-35	ZB-226	A-36
1700-001	A-15	8054-018	A-43	ZB-72	A-37	ZB-144	A-35	ZB-228	A-36
1700-002	A-15	8054-019	A-43	ZB-73	A-37	ZB-145	A-35	ZB-230	A-36
1700-003	A-15	8054-020	A-43	ZB-74	A-37	ZB-147	A-35	ZB-232	A-36
1700-007	A-15	8054-021	A-43	ZB-75	A-37	ZB-148	A-35	ZB-236	A-36
1700-008	A-15	8054-022	A-43	ZB-76	A-37	ZB-149	A-35	ZB-238	A-36
1700-009	A-15	8054-023	A-43	ZB-77	A-37	ZB-152	A-35	ZB-252	A-36
1700-010	A-15	8054-024	A-43	ZB-78	A-37	ZB-154	A-35	ZB-254	A-36
1700-012	A-15	MO-121	A-12	ZB-100	A-34	ZB-160	A-35	ZB-260	A-36
1700-013	A-15	MO-122	A-12	ZB-101	A-34	ZB-162	A-35	ZB-262	A-36
8051-011	A-44	MO-123	A-12	ZB-102	A-34	ZB-164	A-35	ZB-264	A-36
8051-013	A-44	MO-124	A-13	ZB-103	A-34	ZB-168	A-35	ZB-266	A-36
8051-015	A-44	MO-125	A-13	ZB-104	A-34	ZB-170	A-35	ZB-268	A-36
8051-018	A-44	MO-126	A-13	ZB-106	A-34	ZB-180	A-35	ZB-270	A-36
8051-020	A-44	MO-127	A-13	ZB-112	A-34	ZB-181	A-35	ZB-292	A-36
8051-021	A-44	MO-128	A-13	ZB-113	A-34	ZB-182	A-35	ZB-306	A-37
8051-027	A-44	MO-129	A-13	ZB-114	A-34	ZB-183	A-35	ZB-322	A-37
8051-032	A-44	MS-121	A-28	ZB-121	A-34	ZB-184	A-35	ZB-324	A-37

Part #	Page	Part #	Page	Part #	Page	Part #	Page	Part #	Page
ZB-326	A-37	ZB-647	A-38	ZO-067	A-16	ZO-7166	A-7	ZO-7596	A-11
ZB-328	A-37	ZB-648	A-38	ZO-068	A-16	ZO-7167	A-7	ZO-7651	A-11
ZB-330	A-37	ZB-650	A-39	ZO-069	A-16	ZO-7172	A-7	ZO-7654	A-11
ZB-336	A-37	ZB-659	A-39	ZO-070	A-17	ZO-7180	A-8	ZO-7682	A-11
ZB-338	A-37	ZB-662	A-39	ZO-071	A-17	ZO-7201	A-8	ZO-7685	A-11
ZB-360	A-37	ZB-664	A-39	ZO-072	A-17	ZO-7205	A-8	ZO-7689	A-12
ZB-362	A-37	ZO11-089-090	A-3	ZO-073	A-17	ZO-7216	A-8	ZO-7718	A-12
ZB-364	A-37	ZO11-119-090	A-3	ZO-076	A-17	ZO-7225	A-33	ZO-72222	A-3
ZB-366	A-37	ZO11-127-090	A-3	ZO-077	A-17	ZO-7226	A-33	ZO-72622	A-3
ZB-447	A-35	ZO18-089-090	A-3	ZO-079	A-31	ZO-7240	A-8	ZO-72623	A-3
ZB-448	A-35	ZO18-127-090	A-3	ZO-165	A-17	ZO-7241	A-8	ZO-72625	A-3
ZB-449	A-35	ZO18-127-100	A-3	ZO-640	A-5	ZO-7242	A-8	ZO-72689	A-4
ZB-582	A-38	ZO18-135-090	A-3	ZO-7050	A-5	ZO-7243	A-9	ZO-72709	A-4
ZB-584	A-38	ZO25-089-075	A-3	ZO-7052	A-5	ZO-7244	A-9	ZO-72719	A-4
ZB-585	A-38	ZO25-089-090	A-3	ZO-7053	A-6	ZO-7245	A-9	ZO-72730	A-4
ZB-586	A-38	ZO25-107-090	A-3	ZO-7055	A-6	ZO-7246	A-9	ZO-72737	A-4
ZB-617	A-38	ZO25-119-090	A-3	ZO-7079	A-32	ZO-7250	A-9	ZO-73416	A-4
ZB-618	A-38	ZO25-127-090	A-3	ZO-7111	A-6	ZO-7251	A-9	ZO-73417	A-4
ZB-619	A-38	ZO25-127-100	A-3	ZO-7123	A-6	ZO-7274	A-9	ZO-73418	A-4
ZB-625	A-38	ZO25-135-090	A-3	ZO-7125	A-32	ZO-7277	A-10	ZO-73419	A-4
ZB-626	A-38	ZO25-147-090	A-3	ZO-7137	A-32	ZO-7278	A-10	ZO-73420	A-5
ZB-627	A-38	ZO-060	A-16	ZO-7145	A-6	ZO-7302	A-10	ZO-73421	A-5
ZB-628	A-38	ZO-061	A-16	ZO-7151	A-6	ZO-7305	A-10	ZO-73422	A-5
ZB-640	A-39	ZO-062	A-16	ZO-7154	A-7	ZO-7309	A-10	ZO-73489	A-5
ZB-642	A-39	ZO-063	A-16	ZO-7156	A-7	ZO-7377	A-10	ZO-73613	A-5
ZB-644	A-39	ZO-064	A-16	ZO-7161	A-7	ZO-7380	A-10	ZO-73616	A-5
ZB-645	A-39	ZO-065	A-16	ZO-7163	A-7	ZO-7390	A-11	ZR-020	A-30
ZB-646	A-39	ZO-066	A-16	ZO-7165	A-7	ZO-7595	A-11	ZR-021	A-30

Part #	Page	Part #	Page	Part #	Page	Part #	Page
ZR-031	A-31	ZS-029	A-20	ZS-327	A-23	ZS-373	A-19
ZR-032	A-31	ZS-030	A-20	ZS-329	A-23	ZS-382	A-26
ZR-032M	A-31	ZS-031	A-20	ZS-330	A-23	ZS-390	A-26
ZR-033M	A-31	ZS-032	A-21	ZS-331	A-23	ZS-395	A-26
ZR-050	A-30	ZS-033	A-21	ZS-332	A-23	ZS-410	A-25
ZR-051	A-30	ZS-034	A-21	ZS-333	A-23	ZS-412	A-26
ZR-052	A-30	ZS-035	A-21	ZS-334	A-24	ZS-414	A-26
ZR-053	A-30	ZS-036	A-21	ZS-336	A-24	ZS-510	A-26
ZR-058	A-14	ZS-037	A-21	ZS-337	A-24		
ZR-059	A-14	ZS-038	A-21	ZS-338	A-24		
ZR-060	A-14	ZS-038A	A-21	ZS-339	A-24		
ZR-061	A-14	ZS-039	A-22	ZS-340	A-24		
ZR-073	A-14	ZS-040	A-22	ZS-341	A-24		
ZR-079	A-14	ZS-041	A-22	ZS-342	A-24		
ZR-080	A-40	ZS-042	A-22	ZS-343	A-25		
ZR-081	A-40	ZS-043	A-22	ZS-360	A-18		
ZR-082	A-40	ZS-103	A-25	ZS-361	A-18		
ZR-083	A-40	ZS-104	A-22	ZS-362	A-18		
ZR-084	A-40	ZS-105	A-25	ZS-363	A-18		
ZR-085	A-40	ZS-110	A-26	ZS-364	A-18		
ZR-120	A-30	ZS-114	A-25	ZS-365	A-18		
ZR-154	A-30	ZS-115	A-25	ZS-366	A-18		
ZR-160	A-14	ZS-125	A-25	ZS-367	A-18		
ZR-733	A-30	ZS-128	A-22	ZS-368	A-18		
ZS-015	A-20	ZS-129	A-22	ZS-369	A-18		
ZS-018	A-20	ZS-225	A-25	ZS-370	A-19		
ZS-027	A-20	ZS-304	A-23	ZS-371	A-19		
ZS-028	A-20	ZS-318	A-23	ZS-372	A-19		

BRASSELER-STYLE BLADES, BURS + RASPS

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
03-C0312	ZB-154	A-35	04-C0609	ZB-102	A-34	09-D0132	ZB-264	A-36
03-C0901	ZB-136	A-34	04-C0612	ZB-113	A-34	09-D0135	ZB-266	A-36
03-C0905	ZB-138	A-34	04-C0612	ZB-75	A-37	09-D0138	ZB-270	A-36
03-D0135	ZB-164	A-35	04-C0615	ZB-74	A-37	09-D0139	ZB-262	A-36
03-D0138	ZB-168	A-35	04-C0617	ZB-106	A-34	09-D0141	ZB-268	A-36
03-D0911	ZB-170	A-35	04-C1803	ZB-190	A-35	09-D0911	ZB-306	A-37
03-D1901	ZB-180	A-35	04-D0139	ZB-160	A-35	10-C0101	ZB-222	A-36
03-D1902	ZB-181	A-35	04-D1903	ZB-182	A-35	10-C0103	ZB-224	A-36
04-C0101	ZB-121	A-34	04-D2001	ZB-183	A-35	10-C0105	ZB-226	A-36
04-C0103	ZB-122	A-34	04-S1501	ZB-192	A-35	10-C0107	ZB-228	A-36
04-C0104	ZB-76	A-37	04-S1701	ZB-147	A-35	10-C0108	ZB-230	A-36
04-C0105	ZB-124	A-34	04-S1703	ZB-148	A-35	10-C0306	ZB-206	A-36
04-C0106	ZB-71	A-37	04-S1704	ZB-149	A-35	10-C0310	ZB-254	A-36
04-C0107	ZB-126	A-34	06-C0619	ZB-114	A-34	10-C0617	ZB-214	A-36
04-C0108	ZB-128	A-34	06-S1601	ZB-142	A-35	10-D0139	ZB-362	A-37
04-C0109	ZB-130	A-34	07-S1602	ZB-143	A-35	11-C0619	ZB-220	A-36
04-C0170	ZB-132	A-34	07-S1604	ZB-145	A-35	12-C0105	ZB-326	A-37
04-C0302	ZB-101	A-34	07-S1605	ZB-144	A-35	12-C0108	ZB-330	A-37
04-C0303	ZB-78	A-37	08-S1802	ZB-112	A-34	12-C0109	ZB-336	A-37
04-C0306	ZB-104	A-34	09-C0109	ZB-232	A-36	12-C0901	ZB-338	A-37
04-C0307	ZB-72	A-37	09-C0110	ZB-236	A-36	12-C0905	ZB-360	A-37
04-C0310	ZB-152	A-35	09-C0170	ZB-232	A-36	12-D0133	ZB-364	A-37
04-C0404	ZB-77	A-37	09-C0312	ZB-260	A-36	12-D0141	ZB-70	A-37
04-C0602	ZB-100	A-34	09-C0613	ZB-204	A-36	12-S0101	ZB-322	A-37
04-C0604	ZB-103	A-34	09-C0901	ZB-238	A-36	12-S0150	ZB-324	A-37
04-C0606	ZB-73	A-37	09-C0905	ZB-252	A-36	12-S0153	ZB-328	A-37

BRASSELER-STYLE BLADES, BURS + RASPS

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
14-C0104	K-75	B-40*	16-S0109	K-56	B-39*	19-S1301	K-23	B-38*
14-C0106	K-70	B-40*	16-S0110	K-57	B-39*	19-S1705	K-19	B-37*
14-C0308	K-71	B-40*	16-S0111	K-60	B-40*	19-S1802	K-22	B-38*
14-C0606	K-72	B-40*	16-S0113	K-62	B-40*	70-C0108	ZB-619	A-38
14-C0608	K-73	B-40*	16-S0115	K-10	B-37*	70-C0109	ZB-618	A-38
14-C0616	K-74	B-40*	16-S0116	K-64	B-40*	70-C0170	ZB-617	A-38
14-C0702	K-76	B-40*	16-S0117	K-11	B-37*	71-C0107	ZB-628	A-38
14-S0124	K-14	B-37*	16-S0118	K-65	B-40*	71-C0108	ZB-627	A-38
14-S0125	K-15	B-37*	16-S0119	K-52	B-39*	71-C0109	ZB-626	A-38
14-S0127	K-05	B-37*	16-S0120	K-12	B-37*	71-C0170	ZB-625	A-38
14-S0402	K-01	B-37*	16-S0121	K-54	B-39*	71-D0135	ZB-648	A-38
14-S0903	K-02	B-37*	16-S0122	K-55	B-39*	71-D0141	ZB-647	A-38
14-S0909	K-03	B-37*	16-S0123	K-13	B-37*	72-C0107	ZB-646	A-39
14-S1201	K-04	B-37*	16-S0124	K-36	B-38*	72-C0108	ZB-645	A-39
15-S1101	K-08	B-37*	16-S0126	K-58	B-39*	72-C0109	ZB-644	A-39
16-D0128	K-40	B-39*	16-S0402	K-33	B-38*	72-D0133	ZB-659	A-39
16-D0129	K-41	B-39*	16-S0621	K-20	B-38*	72-D0141	ZB-650	A-39
16-D0130	K-42	B-39*	16-S0903	K-38	B-38*	75-D0139	ZB-662	A-39
16-D0131	K-43	B-39*	16-S0909	K-34	B-38*	75-S0101	ZB-664	A-39
16-D0132	K-44	B-39*	16-S1201	K-35	B-38*	75-S0150	ZB-640	A-39
16-D0134	K-45	B-39*	16-S1601	K-26	B-38*	75-S0151	ZB-642	A-39
16-D0135	K-46	B-39*	16-S1702	K-17	B-37*	B60-1190-35	SO11-089-090	B-13*
16-D0136	K-47	B-39*	16-S1801	K-21	B-38*	B60-1190-50	SO11-127-090	B-13*
16-D0137	K-48	B-39*	17-S1602	K-27	B-38*	B60-1890-35	SO18-089-090	B-13*
16-S0103	K-09	B-37*	17-S1603	K-28	B-38*	B60-1890-50	SO18-127-090	B-13*
16-S0106	K-50	B-39*	18-S1701	K-16	B-37*	B60-1890-53	SO18-135-090	B-13*
16-S0107	K-51	B-39*	18-S1703	K-18	B-37*	B60-25100-50	SO25-127-100	B-13*
16-S0108	K-53	B-39*	19-S0801	K-25	B-38*	B60-2575-35	SO25-089-075	B-14*

* Download catalog section B to view product page at www.microaire.com

BRASSELER-STYLE BLADES, BURS + RASPS

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
B60-2590-35	SO25-089-090	B-14*	KM-104	ZS-104	A-22	KM-17002	1700-002	A-15
B60-2590-50	SO25-119-090	B-14*	KM-105	ZS-015	A-20	KM-17003	1700-003	A-15
B60-2590-53	SO25-135-090	B-14*	KM-106	ZS-032	A-21	KM-17007	1700-007	A-15
K121	MS-121	A-28	KM-107	ZS-041	A-22	KM-17008	1700-008	A-15
K122	MS-122	A-28	KM-108	ZS-042	A-22	KM-17009	1700-009	A-15
K123	MS-123	A-28	KM-109	ZS-043	A-22	KM-17010	1700-010	A-15
K124	MS-124	A-29	KM-109C	ZS-030	A-20	KM-17012	1700-012	A-15
K125	MS-125	A-29	KM-12001	1200-001	A-27	KM-17013	1700-013	A-15
K131	MS-131	A-29	KM-12005	1200-005	A-27	KM-218	ZS-018	A-20
K132	MS-132	A-29	KM-12006	1200-006	A-27	KM-21L	MO-121	A-12
K134	MS-134	A-29	KM-12007	1200-007	A-27	KM22L	MO-122	A-12
K-2000-61	SO-677	B-20*	KM-12008	1200-008	A-27	KM-231	ZS-031	A-20
K-2000-65J	SO-699	B-21*	KM-12009	1200-009	A-27	KM-234	ZS-034	A-21
K-2000-65T	SO-676	B-20*	KM-12010	1200-010	A-27	KM-235	ZS-035	A-21
K-2000-69	SO-683	B-20*	KM-12011	1200-011	A-27	KM-237	ZS-037	A-21
KM-000	ZO-062	A-16	KM-12012	1200-012	A-27	KM-23L	MO-123	A-12
KM-001	ZO-064	A-16	KM-12013	1200-013	A-28	KM-240	ZS-040	A-22
KM-002	ZO-073	A-17	KM-12014	1200-014	A-28	KM-24L	MO-124	A-13
KM-004	ZO-065	A-16	KM-12017	1200-017	A-28	KM-250R	SR-250	B-22*
KM-005	ZO-063	A-16	KM-12018	1200-018	A-28	KM-25L	MO-125	A-13
KM-006	ZO-060	A-16	KM-12019	1200-019	A-28	KM-261	ZO-061	A-16
KM-008	ZO-066	A-16	KM-14001	1400-001	A-29	KM-267	ZO-067	A-16
KM-060	ZR-020	A-30	KM-14002	1400-002	A-29	KM-268	ZO-068	A-16
KM-070	ZR-021	A-30	KM-14003	1400-003	A-29	KM-269	ZO-069	A-16
KM-100C	ZS-033	A-31	KM-14004	1400-004	A-29	KM-26L	MO-126	A-13
KM-101	ZS-038	A-21	KM-14005	1400-005	A-29	KM-275R	SR-275	B-22*
KM-101C	ZS-038A	A-21	KM-14006	1400-006	A-29	KM-277	SR-027	B-36*
KM-102	ZS-039	A-22	KM-17001	1700-001	A-15	KM-279	ZO-079	A-31

* Download catalog section B to view product page at www.microaire.com

BRASSELER-STYLE BLADES, BURS + RASPS

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
KM-27L	MO-127	A-13	KM3-115	SP-115	B-26*	KM33-125	SP-125A	B-27*
KM-28L	MO-128	A-13	KM3-125	SP-125	B-27*	KM33-206	SP-206A	B-28*
KM-29L	MO-129	A-13	KM-3170	ZS-370	A-19	KM33-212	SP-212A	B-28*
KM-3000	ZS-362	A-18	KM-3171	ZS-371	A-19	KM33-225	SP-225A	B-28*
KM-3001	ZS-364	A-18	KM-3172	ZS-372	A-19	KM33-410	SP-410A	B-29*
KM-3002	ZS-369	A-18	KM-3173	ZS-373	A-19	KM33-411	SP-411A	B-30*
KM-3005	ZS-363	A-18	KM-32	ZR-032	A-31	KM33-412	SP-412A	B-30*
KM-3006	ZS-360	A-18	KM3-206	SP-206	B-27*	KM33-413	SP-413A	B-30*
KM-3008	ZS-366	A-18	KM3-212	SP-212	B-28*	KM33-414	SP-414A	B-30*
KM-300R	SR-300	B-22*	KM-3218	ZS-318	A-23	KM3-410	SP-410	B-29*
KM31	ZR-031	A-31	KM3-225	SP-225	B-28*	KM3-411	SP-411	B-29*
KM-3100C	ZS-333	A-23	KM-3229	ZS-329	A-23	KM3-412	SP-412	B-30*
KM-3101	ZS-338	A-24	KM-3231	ZS-331	A-23	KM3-413	SP-413	B-30*
KM3-101	SP-101	B-25*	KM-3237	ZS-337	A-24	KM3-414	SP-414	B-30*
KM-3101C	ZS-343	A-25	KM-3240	ZS-340	A-24	KM410S	SO-410	B-23*
KM-3102	ZS-339	A-24	KM-325R	SR-325	B-22*	KM-411S	SO-411	B-23*
KM3-102	SP-102	B-25*	KM-3261	ZS-361	A-18	KM-412S	SO-412	B-23*
KM3-103	SP-103	B-25*	KM-3265	ZS-365	A-18	KM-413S	SO-413	B-23*
KM3-104	SP-104	B-25*	KM33-101	SP-101A	B-25*	KM-414S	SO-414	B-23*
KM3-105	SP-105	B-25*	KM33-102	SP-102A	B-25*	KM-415S	SO-415	B-24*
KM-3106	ZS-332	A-23	KM33-103	SP-103A	B-25*	KM-416S	SO-416	B-24*
KM3-106	SP-106	B-26*	KM33-104	SP-104A	B-25*	KM-417S	SO-417	B-24*
KM-3107	ZS-341	A-24	KM33-105	SP-105A	B-26*	KM-458	ZR-058	A-14
KM-3108	ZS-342	A-24	KM33-106	SP-106A	B-26*	KM-459	ZR-059	A-14
KM3-108	SP-108	B-26*	KM33-108	SP-108A	B-26*	KM-460	ZR-060	A-14
KM-3109C	ZS-330	A-23	KM33-111	SP-111A	B-26*	KM-461	ZR-061	A-14
KM3-111	SP-111	B-26*	KM33-114	SP-114A	B-26*	KM-473	ZR-073	A-14
KM3-114	SP-114	B-26*	KM33-115	SP-115A	B-27*	KM-479	ZR-079	A-14

* Download catalog section B to view product page at www.microaire.com

BRASSELER-STYLE BLADES, BURS + RASPS

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
KM-50	ZR-050	A-30	KM-656	SO-656	B-19*	KM71-586	ZO-7165	A-7
KM-51	ZR-051	A-30	KM-658	SO-658	B-19*	KM71-588	ZO-7278	A-10
KM5100-137-115	SR-51-080	B-42*	KM-660	SO-660	B-19*	KM71-595	ZO-7595	A-11
KM5100-137-123	SR-51-723	B-33*	KM-675	SO-675	B-19*	KM71-596	ZO-7596	A-11
KM5100-137-133	SR-51-733	B-33*	KM-683	SO-682	B-20*	KM-72	ZO-072	A-17
KM5100-137-233	SR-51-007	B-32*	KM-689	SO-685	B-20*	KM-7201	ZO-7201	A-8
KM5100-37-113	SR-51-082	B-42*	KM-693	SO-689	B-20*	KM-7567	ZO-7167	A-7
KM5100-37-114	SR-51-083	B-42*	KM-70	ZO-070	A-17	KM-7572	ZO-7172	A-7
KM5100-37-116	SR-51-081	B-42*	KM-71	ZO-071	A-17	KM-80	ZR-080	A-40
KM-52	ZR-052	A-30	KM71-523	ZO-7123	A-6	KM-801	SR-080	B-41*
KM-53	ZR-053	A-30	KM71-540	ZO-7240	A-8	KM-802	SR-083	B-41*
KM-600	SO-600	B-15*	KM71-542	ZO-7242	A-8	KM-81	ZR-081	A-40
KM-602	SO-602	B-15*	KM71-543	ZO-7243	A-9	KM-82	ZR-082	A-40
KM-605	SO-605	B-15*	KM71-544	ZO-7244	A-9	KM-83	ZR-083	A-40
KM-609	SO-609	B-15*	KM71-545	ZO-7245	A-9	KM-97	SR-100	B-36*
KM-610	SO-610	B-16*	KM71-550	ZO-7050	A-9	KM-98	SR-101	B-36*
KM-611	SO-611	B-16*	KM71-551	ZO-7151	A-6	ST2108-100	SO6-100	B-3*
KM-615	SO-615	B-16*	KM71-552	ZO-7052	A-5	ST2108-102	SO6-102	B-3*
KM-618	SO-618	B-16*	KM71-553	ZO-7053	A-6	ST2108-103	SO6-103	B-3*
KM-619	SO-718	B-21*	KM71-554	ZO-7154	A-7	ST2108-105	SO6-105	B-3*
KM-620	SO-620	B-16*	KM71-555	ZO-7055	A-6	ST2108-109	SO6-109	B-3*
KM-625	SO-625	B-35*	KM71-556	ZO-7156	A-7	ST2108-110	SO6-110	B-4*
KM-640	SO-640	B-18*	KM71-561	ZO-7161	A-7	ST2108-111	SO6-111	B-4*
KM-645	SO-645	B-18*	KM71-563	ZO-7163	A-7	ST2108-113	SO6-113	B-4*
KM-648	SO-648	B-18*	KM71-566	ZO-7166	A-7	ST2108-115	SO6-115	B-4*
KM-650	SO-650	B-18*	KM71-574	ZO-7274	A-9	ST2108-118	SO6-118	B-4*
KM-651	SO-651	B-18*	KM71-577	ZO-7277	A-10	ST2108-120	SO6-120	B-4*
KM-655	SO-655	B-19*	KM71-580	ZO-7180	A-8	ST2108-125	SO6-125	B-34*

* Download catalog section B to view product page at www.microaire.com

BRASSELER-STYLE BLADES, BURS + RASPS

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
ST2108-131	SO6-131	B-4*	ST4118-127-100	SO618-127-100	B-11*
ST2108-140	SO6-140	B-4*	ST4118-137-090	SO618-135-090	B-11*
ST2108-144	SO6-144	B-5*	ST4125-064-075	SO625-064-075	B-11*
ST2108-145	SO6-145	B-5*	ST4125-089-090	SO625-089-075	B-11*
ST2108-148	SO6-148	B-5*	ST4125-127-090	SO625-127-090	B-12*
ST2108-150	SO6-150	B-5*	ST4125-127-100	SO625-127-100	B-12*
ST2108-151	SO6-151	B-5*	ST4125-135-090	SO625-135-090	B-12*
ST2108-152	SO6-152	B-5*	ST4125-147-090	SO625-147-090	B-12*
ST2108-155	SO6-155	B-5*			
ST2108-156	SO6-156	B-6*			
ST2108-158	SO6-158	B-6*			
ST2108-160	SO6-160	B-6*			
ST2108-175	SO6-175	B-6*			
ST2108-176	SO6-176	B-6*			
ST2108-177	SO6-177	B-6*			
ST2108-182	SO6-182	B-7*			
ST2108-183	SO6-183	B-7*			
ST2108-185	SO6-185	B-7*			
ST2108-189	SO6-189	B-7*			
ST2108-193	SO6-193	B-7*			
ST2108-195	SO6-195	B-8*			
ST2108-197	SO6-197	B-8*			
ST2108-218	SO6-218	B-8*			
ST4111-127-090	SO611-127-090	B-10*			
ST4111-147-090	SO611-147-075	B-10*			
ST4118-089-075	SO618-089-075	B-11*			
ST4118-119-090	SO618-119-090	B-11*			
ST4118-127-090	SO618-127-075	B-11*			

* Download catalog section B to view product page at www.microaire.com

BRASSELER-STYLE DRILLS, PINS + WIRES

Brasseler	MicroAire	Page	Brasseler	MicroAire	Page	Brasseler	MicroAire	Page
KM166-00-00	8054-010NS	A-43	KM168-19-31	1648-509NS	A-42	KM169-19-96	1636-509TNS	A-42
KM166-00-01	8054-011NS	A-43	KM168-19-33	1624-509NS	A-42	KM169-29-18	1632-109TNS	A-42
KM166-00-02	8054-012NS	A-43	KM168-19-53	1640-509NS	A-42	KM169-29-31	1648-109TNS	A-42
KM166-00-03	8054-013NS	A-43	KM168-19-56	1620-509NS	A-42	KM169-29-33	1624-109TNS	A-42
KM166-00-04	8054-014NS	A-43	KM168-19-76	1628-509NS	A-42	KM169-29-53	1640-109TNS	A-42
KM166-00-05	8054-015NS	A-43	KM168-19-96	1636-509NS	A-42	KM169-29-56	1620-109TNS	A-42
KM166-00-06	8054-016NS	A-43	KM168-29-18	1632-109NS	A-42	KM169-29-76	1628-109TNS	A-42
KM166-00-07	8054-018NS	A-43	KM168-29-31	1648-109NS	A-42	KM169-29-96	1636-109TNS	A-42
KM166-00-08	8054-019NS	A-43	KM168-29-33	1624-109NS	A-42	KM169-39-18	1632-609TNS	A-42
KM166-310-11	8053-011NS	C-4*	KM168-29-53	1640-109NS	A-42	KM169-39-31	1648-609TNS	A-42
KM166-310-15	8053-015NS	C-4*	KM168-29-56	1620-109NS	A-42	KM169-39-33	1624-609TNS	A-42
KM166-310-16	8053-016NS	C-4*	KM168-29-76	1628-109NS	A-42	KM169-39-53	1640-609TNS	A-42
KM166-310-19	8053-019NS	C-4*	KM168-29-96	1636-109NS	A-42	KM-169-39-56	1620-609TNS	A-42
KM166-310-21	8053-021NS	C-4*	KM168-39-18	1632-609NS	A-42	KM169-39-76	1628-609TNS	A-42
KM166-310-23	8053-023NS	C-4*	KM168-39-31	1648-609NS	A-42	KM169-39-96	1636-609TNS	A-42
KM166-310-25	8053-025NS	C-4*	KM168-39-33	1624-609NS	A-42	KM-172-155-28	1600-020NS	A-41
KM166-310-26	8053-026NS	C-4*	KM168-39-53	1640-609NS	A-42	KM-172-155-35	1600-021NS	A-41
KM166-310-28	8053-028NS	C-4*	KM168-39-56	1620-609NS	A-42	KM-172-155-45	1600-022NS	A-41
KM166-310-29	8053-029NS	C-4*	KM168-39-76	1628-609NS	A-42	KM-172-155-62	1600-023NS	A-41
KM166-310-31	8053-031NS	C-4*	KM168-39-96	1636-609NS	A-42	KM-172-19-35	1600-9355NS	A-41
KM166-310-35	8053-035NS	C-4*	KM169-19-18	1632-509TNS	A-42	KM-172-19-45	1600-9455NS	A-41
KM166-310-37	8053-037NS	C-4*	KM169-19-31	1648-509TNS	A-42	KM-172-19-54	1600-9545NS	A-41
KM166-310-44	8053-044NS	C-4*	KM169-19-33	1624-509TNS	A-42	KM-172-19-62	1600-9625NS	A-41
KM166-310-48	8053-048NS	C-4*	KM169-19-53	1640-509TNS	A-42	KM-172-24-28	1600-428NS	A-41
KM166-310-60	8053-060NS	C-4*	KM169-19-56	1620-509TNS	A-42	KM-172-24-35	1600-435NS	A-41
KM168-19-18	1632-509NS	A-42	KM169-19-76	1628-509TNS	A-42	KM-172-24-45	1600-445NS	A-41

* Download catalog section C to view product page at www.microaire.com

BRASSELER-STYLE DRILLS, PINS + WIRES

Brasseler	MicroAire	Page
KM-172-24-54	1600-454NS	A-41
KM-172-24-62	1600-462NS	A-41
KM-172-26-28	1600-628NS	A-41
KM-172-26-35	1600-635NS	A-41
KM-172-26-45	1600-645NS	A-41
KM-172-26-54	1600-654NS	A-41
KM-172-26-62	1600-662NS	A-41
KM-172-29-35	1600-935NS	A-41
KM-172-29-45	1600-945NS	A-41
KM-172-29-54	1600-954NS	A-41
KM-172-29-62	1600-962NS	A-41
KM-173-29-35	1600-935TNS	A-41
KM-173-29-45	1600-945TNS	A-41
KM-173-29-54	1600-954TNS	A-41
KM-173-29-62	1600-962TNS	A-41

LINVATEC-STYLE BLADES, BURS + RASPS

Linvatec	MicroAire	Page	Linvatec	MicroAire	Page	Linvatec	MicroAire	Page
1384-77	ZO-077	A-17	5023-171	ZS-371	A-19	5053-037/237	ZS-037	A-21
5021-163	ZS-363	A-18	5023-172	ZS-372	A-19	5053-038/238	ZS-038	A-21
5023-118	ZS-318	A-23	5023-173	ZS-373	A-19	5053-039/239	ZS-039	A-22
5023-129	ZS-329	A-23	5023-179	ZO-079	A-31	5053-040/240	ZS-040	A-22
5023-130	ZS-330	A-23	5052-173	ZR-073	A-14	5053-041/241	ZS-041	A-22
5023-131	ZS-331	A-23	5052-179	ZR-079	A-14	5053-042/242	ZS-042	A-22
5023-132	ZS-332	A-23	5052-182	ZB-582	A-38	5053-050	ZR-050	A-30
5023-133	ZS-333	A-23	5052-184	ZB-584	A-38	5053-051	ZR-051	A-30
5023-134	ZS-334	A-24	5052-185	ZB-585	A-38	5053-052	ZR-052	A-30
5023-136	ZS-336	A-24	5052-186	ZB-586	A-38	5053-053	ZR-053	A-30
5023-137	ZS-337	A-24	5052-258	ZR-058	A-14	5053-054	ZR-154	A-30
5023-138	ZS-338	A-24	5052-259	ZR-059	A-14	5053-060/260	ZO-060	A-16
5023-139	ZS-339	A-24	5052-260	ZR-160	A-14	5053-061/261	ZO-061	A-16
5023-140	ZS-340	A-24	5052-261	ZR-060	A-14	5053-062/262	ZO-062	A-16
5023-141	ZS-341	A-24	5052-261	ZR-061	A-14	5053-063/263	ZO-063	A-16
5023-142	ZS-342	A-24	5052-279	ZR-079	A-14	5053-064/264	ZO-064	A-16
5023-143	ZS-343	A-25	5053-018/218	ZS-018	A-20	5053-067/267	ZO-067	A-16
5023-160	ZS-360	A-18	5053-028/228	ZS-128	A-22	5053-069/269	ZO-069	A-16
5023-161	ZS-361	A-18	5053-029/229	ZS-129	A-22	5053-073/273	ZO-073	A-17
5023-162	ZS-362	A-18	5053-030/230	ZS-030	A-20	5053-170	ZO-070	A-17
5023-163	ZS-363	A-18	5053-031/231	ZS-031	A-20	5053-171	ZO-071	A-17
5023-164	ZS-364	A-18	5053-032/232	ZS-032	A-21	5053-172	ZO-072	A-17
5023-165	ZS-365	A-18	5053-033/233	ZS-033	A-21	5053-220	ZR-020	A-30
5023-166	ZS-366	A-18	5053-034/234	ZS-034	A-21	5053-220	ZR-120	A-30
5023-169	ZS-369	A-18	5053-035/235	ZS-035	A-21	5053-221	ZR-021	A-30
5023-170	ZS-370	A-19	5053-036/236	ZS-036	A-21	5053-80	ZR-080	A-40

LINVATEC-STYLE BLADES, BURS + RASPS

Linvatec	MicroAire	Page	Linvatec	MicroAire	Page	Linvatec	MicroAire	Page
5053-81	ZR-081	A-40	5071-140	ZO-7240	A-8	5071-572	ZO-7172	A-7
5053-82	ZR-082	A-40	5071-141	ZO-7241	A-8	5091-080	ZB-180	A-35
5053-83	ZR-083	A-40	5071-142	ZO-7242	A-8	5091-081	ZB-181	A-35
5053-84	ZR-084	A-40	5071-143	ZO-7243	A-9	5091-082	ZB-182	A-35
5056-117	ZB-617	A-38	5071-144	ZO-7244	A-9	5091-083	ZB-183	A-35
5056-118	ZB-618	A-38	5071-145	ZO-7245	A-9	5091-201	ZB-101	A-34
5056-119	ZB-619	A-38	5071-146	ZO-7246	A-9	5091-202	ZB-102	A-34
5056-125	ZB-625	A-38	5071-161	ZO-7161	A-7	5091-204	ZB-104	A-34
5056-126	ZB-626	A-38	5071-163	ZO-7163	A-7	5091-206	ZB-106	A-34
5056-127	ZB-627	A-38	5071-167	ZO-7167	A-7	5091-212	ZB-112	A-34
5056-128	ZB-628	A-38	5071-174	ZO-7274	A-9	5091-214	ZB-114	A-34
5056-140	ZB-640	A-39	5071-180	ZO-7180	A-8	5091-221	ZB-121	A-34
5056-142	ZB-642	A-39	5071-181	ZO-7201	A-8	5091-222	ZB-122	A-34
5056-144	ZB-644	A-39	5071-186	ZO-7165	A-7	5091-224	ZB-124	A-34
5056-145	ZB-645	A-39	5071-186	ZO-7166	A-7	5091-225	ZB-125	A-34
5056-146	ZB-646	A-39	5071-187	ZO-7277	A-10	5091-226	ZB-126	A-34
5056-147	ZB-647	A-38	5071-188	ZO-7278	A-10	5091-228	ZB-128	A-34
5056-148	ZB-648	A-38	5071-195	ZO-7595	A-11	5091-230	ZB-130	A-34
5056-150	ZB-650	A-39	5071-196	ZO-7596	A-11	5091-232	ZB-132	A-34
5056-159	ZB-659	A-39	5071-201	ZO-7201	A-8	5091-236	ZB-136	A-34
5056-162	ZB-662	A-39	5071-530	ZO-7390	A-11	5091-238	ZB-138	A-34
5056-164	ZB-664	A-39	5071-550	ZO-7050	A-5	5091-242	ZB-142	A-35
5059-031	ZR-031	A-31	5071-551	ZO-7151	A-6	5091-243	ZB-143	A-35
5059-031	ZR-032	A-31	5071-552	ZO-7052	A-5	5091-244	ZB-144	A-35
5059-279	ZO-7079	A-32	5071-553	ZO-7053	A-6	5091-245	ZB-145	A-35
5071-123	ZO-7123	A-6	5071-554	ZO-7154	A-7	5091-247	ZB-147	A-35
5071-125	ZO-7225	A-33	5071-555	ZO-7055	A-6	5091-248	ZB-148	A-35
5071-126	ZO-7226	A-33	5071-556	ZO-7156	A-7	5091-249	ZB-149	A-35

LINVATEC-STYLE BLADES, BURS + RASPS

Linvatec	MicroAire	Page	Linvatec	MicroAire	Page	Linvatec	MicroAire	Page
5091-252	ZB-152	A-35	5093-266	ZB-366	A-37	KM-17009	1700-009	A-15
5092-203	ZB-203	A-36	K121	MS-121	A-28	KM-17010	1700-010	A-15
5092-206	ZB-206	A-36	K122	MS-122	A-28	KM-17012	1700-012	A-15
5092-214	ZB-214	A-36	K123	MS-123	A-28	KM-17013	1700-013	A-15
5092-220	ZB-220	A-36	K124	MS-124	A-29	L121	MO-121	A-12
5092-222	ZB-222	A-36	K125	MS-125	A-29	L122	MO-122	A-12
5092-224	ZB-224	A-36	K131	MS-131	A-29	L123	MO-123	A-12
5092-226	ZB-226	A-36	K132	MS-132	A-29	L124	MO-124	A-13
5092-228	ZB-228	A-36	K134	MS-134	A-29	L125	MO-125	A-13
5092-230	ZB-230	A-36	KM-12001	1200-001	A-27	L126	MO-126	A-13
5092-232	ZB-232	A-36	KM-12005	1200-005	A-27	L127	MO-127	A-13
5092-236	ZB-236	A-36	KM-12006	1200-006	A-27	L128	MO-128	A-13
5092-238	ZB-238	A-36	KM-12007	1200-007	A-27	L129	MO-129	A-13
5092-252	ZB-252	A-36	KM-12008	1200-008	A-27	ZR-032M	ZR-032M	A-31
5092-254	ZB-154	A-35	KM-12009	1200-009	A-27	ZR-033M	ZR-033M	A-31
5092-254	ZB-254	A-36	KM-12010	1200-010	A-27	ZS-110	ZS-110	A-26
5092-260	ZB-260	A-36	KM-12011	1200-011	A-27	ZS-304	ZS-304	A-23
5092-264	ZB-264	A-36	KM-12012	1200-012	A-27	ZS-327	ZS-327	A-23
5092-268	ZB-268	A-36	KM-12013	1200-013	A-28	ZS-382	ZS-382	A-26
5092-270	ZB-270	A-36	KM-12014	1200-014	A-28	ZS-390	ZS-390	A-26
5093-222	ZB-322	A-37	KM-12017	1200-017	A-28	ZS-395	ZS-395	A-26
5093-224	ZB-324	A-37	KM-12018	1200-018	A-28	ZS-510	ZS-510	A-26
5093-226	ZB-326	A-37	KM-12019	1200-019	A-28			
5093-228	ZB-328	A-37	KM-17001	1700-001	A-15			
5093-230	ZB-330	A-37	KM-17002	1700-002	A-15			
5093-236	ZB-336	A-37	KM-17003	1700-003	A-15			
5093-238	ZB-338	A-37	KM-17007	1700-007	A-15			
5093-262	ZB-362	A-37	KM-17008	1700-008	A-15			

LINVATEC-STYLE DRILLS, PINS + WIRES

Linvatec	MicroAire	Page	Linvatec	MicroAire	Page	Linvatec	MicroAire	Page
186-04-16	1600-628NS	A-41	187-05-69	1636-609NS	A-42	2410-45	8053-044	C-4*
186-44	1600-635NS	A-41	187-05-79	1636-709NS	A-42	2410-46	8053-048	C-4*
186-46	1600-935NS	A-41	187-06-19	1640-109NS	A-42	2410-600	8053-060	C-4*
186-60	1600-645NS	A-41	187-06-59	1640-509NS	A-42	256-01-21	8054-010	A-43
186-62	1600-945NS	A-41	187-06-69	1640-609NS	A-42	256-01-22	8054-011	A-43
186-66	1600-9455NS	A-41	187-06-79	1640-709NS	A-42	256-01-23	8054-012	A-43
186-76	1600-662NS	A-41	187-07-19	1648-109NS	A-42	256-01-24	8054-013	A-43
186-78	1600-962NS	A-41	187-07-59	1648-509NS	A-42	256-01-25	8054-014	A-43
186-82	1600-9625NS	A-41	187-07-69	1648-609NS	A-42	256-01-26	8054-015	A-43
187-01-19	1620-109NS	A-42	187-07-79	1648-709NS	A-42	256-01-27	8054-016	A-43
187-01-69	1620-609NS	A-42	2318-20	8054-020	A-43	256-01-28	8054-017	A-43
187-01-79	1620-709NS	A-42	2318-27	8054-021	A-43	256-01-29	8054-018	A-43
187-02-19	1624-109NS	A-42	2318-32	8054-022	A-43	2620-01	1636-109TNS	A-42
187-02-59	1624-509NS	A-42	2318-35	8054-023	A-43	2620-03	1628-109TNS	A-42
187-02-69	1624-609NS	A-42	2318-45	8054-024	A-43	262-01-04	1648-609TNS	A-42
187-02-79	1624-709NS	A-42	2410-11	8053-011	C-4*	2620-11	1648-109TNS	A-42
187-03-19	1628-109NS	A-42	2410-15	8053-015	C-4*	2620-13	1648-509TNS	A-42
187-03-59	1628-509NS	A-42	2410-20	8053-019	C-4*	2620-15	1640-109TNS	A-42
187-03-69	1628-609NS	A-42	2410-21	8053-021	C-4*	2620-17	1640-509TNS	A-42
187-03-79	1628-709NS	A-42	2410-25	8053-025	C-4*	262-01-79	1648-709TNS	A-42
187-04-19	1632-109NS	A-42	2410-27	8053-026	C-4*	2620-19	1636-509TNS	A-42
187-04-59	1632-509NS	A-42	2410-28	8053-028	C-4*	262-02-04	1640-609TNS	A-42
187-04-69	1632-609NS	A-42	2410-32	8053-031	C-4*	2620-21	1632-109TNS	A-42
187-04-79	1632-709NS	A-42	2410-33	8053-029	C-4*	2620-23	1632-509TNS	A-42
187-05-19	1636-109NS	A-42	2410-35	8053-035	C-4*	2620-25	1628-509TNS	A-42
187-05-59	1636-509NS	A-42	2410-36	8053-037	C-4*	2620-27	1624-109TNS	A-42

* Download catalog section C to view product page at www.microaire.com

LINVATEC-STYLE DRILLS, PINS + WIRES

Linvatec	MicroAire	Page	Linvatec	MicroAire	Page	Linvatec	MicroAire	Page
262-02-79	1640-709TNS	A-42	47-186-76	1600-662	A-41	47-187-07-59	1648-509	A-42
2620-29	1624-509TNS	A-42	47-186-78	1600-962	A-41	47-187-07-69	1648-609	A-42
262-03-04	1636-609TNS	A-42	47-186-82	1600-9625	A-41	47-187-07-79	1648-709	A-42
2620-31	1620-109TNS	A-42	47-187-01-19	1620-109	A-42	47-2620-01	1636-109T	A-42
2620-33	1620-509TNS	A-42	47-187-01-59	1620-509	A-42	47-2620-03	1628-109T	A-42
2620-35	1600-962TNS	A-41	47-187-01-69	1620-609	A-42	47-262-01-04	1648-609T	A-42
2620-37	1600-945TNS	A-41	47-187-01-79	1620-709	A-42	47-2620-11	1648-109T	A-42
262-03-79	1636-709TNS	A-42	47-187-02-19	1624-109	A-42	47-2620-13	1648-509T	A-42
2620-39	1600-935TNS	A-41	47-187-02-59	1624-509	A-42	47-2620-15	1640-109T	A-42
262-05-04	1628-609TNS	A-42	47-187-02-69	1624-609	A-42	47-2620-17	1640-509T	A-42
262-05-79	1628-709TNS	A-42	47-187-02-79	1624-709	A-42	47-262-01-79	1648-709T	A-42
262-06-04	1624-609TNS	A-42	47-187-03-19	1628-109	A-42	47-2620-19	1636-509T	A-42
262-06-04	1632-609TNS	A-42	47-187-03-59	1628-509	A-42	47-262-02-04	1640-609T	A-42
262-06-79	1624-709TNS	A-42	47-187-03-69	1628-609	A-42	47-2620-21	1632-109T	A-42
262-06-79	1632-709TNS	A-42	47-187-03-79	1628-709	A-42	47-2620-23	1632-509T	A-42
262-07-04	1620-609TNS	A-42	47-187-04-19	1632-109	A-42	47-2620-25	1628-509T	A-42
262-07-79	1620-709TNS	A-42	47-187-04-59	1632-509	A-42	47-2620-27	1624-109T	A-42
286-50	1600-9355NS	A-41	47-187-04-69	1632-609	A-42	47-262-02-79	1640-709T	A-42
287-01-59	1620-509NS	A-42	47-187-04-79	1632-709	A-42	47-2620-29	1624-509T	A-42
47-186-04-14	1600-428	A-41	47-187-05-19	1636-109	A-42	47-262-03-04	1636-609T	A-42
47-186-04-16	1600-628	A-41	47-187-05-59	1636-509	A-42	47-2620-31	1620-109T	A-42
47-186-04-54	1604-128	A-41	47-187-05-69	1636-609	A-42	47-2620-33	1620-509T	A-42
47-186-44	1600-635	A-41	47-187-05-79	1636-709	A-42	47-2620-35	1600-962T	A-41
47-186-46	1600-935	A-41	47-187-06-19	1640-109	A-42	47-2620-37	1600-945T	A-41
47-186-50	1600-9355	A-41	47-187-06-59	1640-509	A-42	47-262-03-79	1636-709T	A-42
47-186-60	1600-645	A-41	47-187-06-69	1640-609	A-42	47-2620-39	1600-935T	A-41
47-186-62	1600-945	A-41	47-187-06-79	1640-709	A-42	47-262-05-04	1628-609T	A-42
47-186-66	1600-9455	A-41	47-187-07-19	1648-109	A-42	47-262-05-79	1628-709T	A-42

LINVATEC-STYLE DRILLS, PINS + WIRES

Linvatec	MicroAire	Page
47-262-06-04	1624-609T	A-42
47-262-06-04	1632-609T	A-42
47-262-06-79	1624-709T	A-42
47-262-06-79	1632-709T	A-42
47-262-07-04	1620-609T	A-42
47-262-07-79	1620-709T	A-42

All dimensions shown are nominal values subject to standard variances and tolerances.

MicroAire® is a registered trademark of MicroAire Surgical Instruments LLC.

The following companies are not affiliated with MicroAire Surgical Instruments LLC:

Conmed Corporation

Dyonics, Inc

Harris/Aesculap

Jacobs

Linvatec Corporation

Stryker

Synthes (USA)

Terumo

Zimmer, Inc

Dyonics® is a trademark of Smith&Nephew

Jacobs® is a registered trademark of Jacobs

Linvatec®, Hall®, Cebotome™, Osteon™, and Neurairtome™ Mini-Driver®, Maxi-Driver® and Maxion™ are trademarks of Linvatec® Corporation, a subsidiary of ConMed Corporation.

Stryker® is a registered trademark of Stryker Corporation

Synthes® is a registered trademark of DePuy Synthes

Terumo® is a registered trademark of Terumo Medical Corporation

Zimmer® is a registered trademark of Zimmer, Inc

MICROAIRE®
For Surgery. For Life.™

MicroAire Surgical Instruments

3590 Grand Forks Boulevard

Charlottesville, VA 22911 USA

800 722 0822 or +434 975 8000

www.microaire.com

© 2016 MicroAire Surgical Instruments, LLC
All rights reserved. LIT-DISPOSABLES-CAT REV A



A COLSON ASSOCIATE