

HealFix Ti Suture Anchors Family

Various options in anchor design, diameter, suture size/color, with/without needles



reddot winner 2020

The Red Dot Award:Product Design only awards this sought-after seal of quality to products that feature an outstanding design. In 2020, HealFix Suture Anchor won the Red Dot award in the category“Medical Devices and Technology”.

-  1.3mm with needles
2.0mm with needles
2.4mm with needles
-  2.0mm with needles
2.0mm
-  2.8mm with needles
2.8mm
-  3.5mm with needles
3.5mm
-  5.0mm with needles
5.0mm
6.5mm with needles
6.5mm
-  4.5mm Plus
5.5mm Plus
6.5mm Plus



Hangzhou Rejoin Mastin Medical Device Co., Ltd

Headquarter Office: 22 Xinyan Rd, Linping Economic Development Zone, Hangzhou, Zhejiang, 311100, P.R. of China
Beijing Branch: Room 1605, Jinghui building, No. 118 Jianguo Road, Chaoyang District, Beijing, P.R. of China
Shanghai Branch: 309A, L1-B, Hongqiao World Center, 1588 Zhuguang Rd, Qingpu District, Shaghai, P.R. of China

contact@rejoin.com @rejoinmedical www.rejoin.com
www.rejoin-medical.com

The content in this document is provided for reference only




HealFix Titanium
Anchor Fixation System

Wide Range of Applications

More options for sizes, threads & sutures

Pre-loaded with WiLSuture Ultra sutures

Ergonomic handle design provides surgeons comfort and ease of use

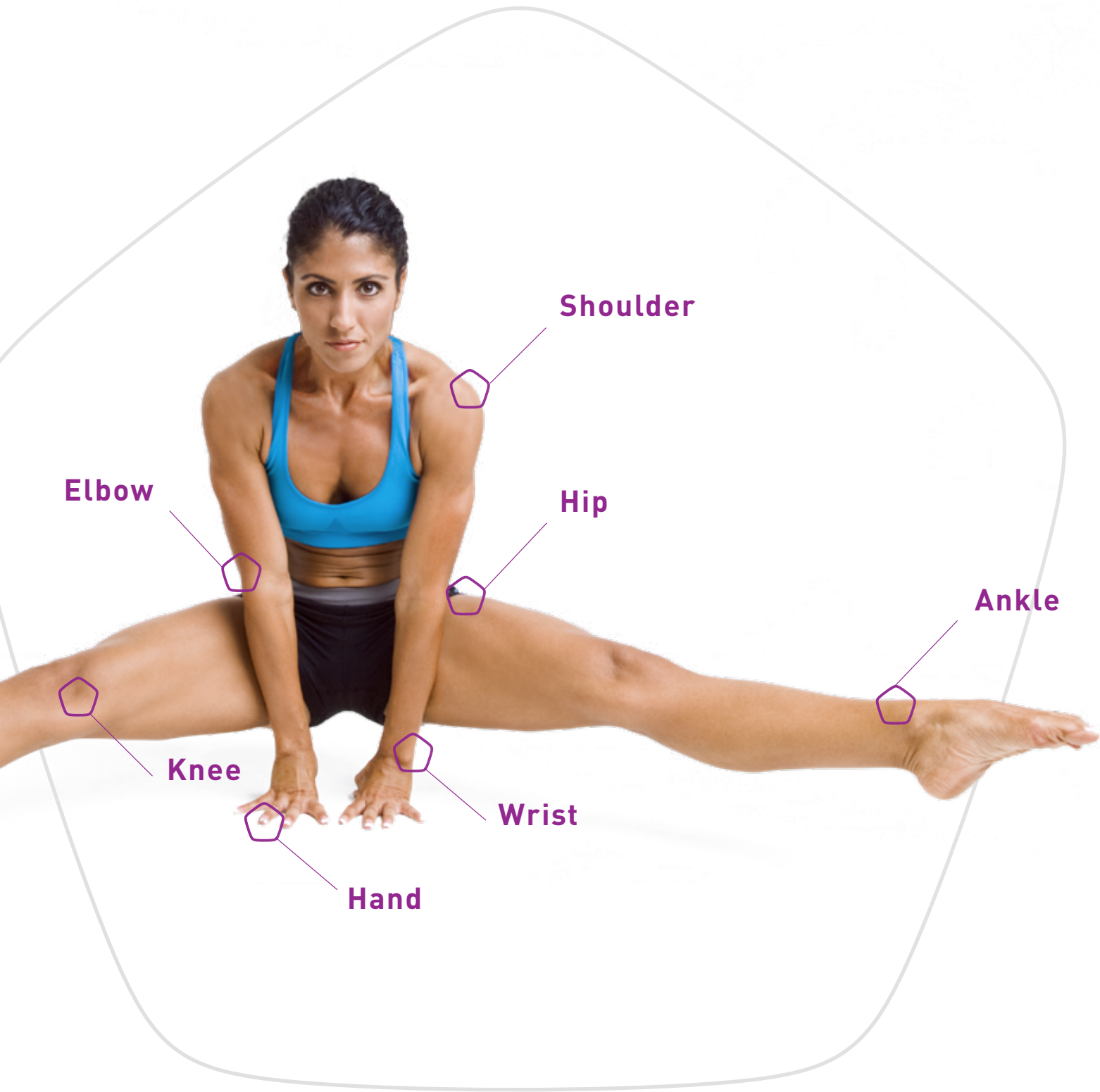

HealFix Ti Suture Anchor

- High-Low double thread design enhances the pull-out strength
- Two independent eyelets provide optimal fixation of both sutures & eliminate suture twining
- Sharpened tip facilitates easy insertion



With needles

Without needles

HealFix Plus Ti Suture Anchor

- Fully Helix transitional thread design increases contact areas with cortical bone for high pull-out strength
- Built-in suture hang-up hole eliminates Cheese-Wire effect between sutures and bone after screw implantation
- Self-punching design makes implantation easier



HealFix Mini Ti Suture Anchor

- NiTi arcs design provides subcortical purchase with outstanding pull-out strength
- Smaller size allows comprehensive fixation and restoration of natural footprint
- 1.3mm minimum size provides optimized soft tissue fixation for DEX applications



Ordering Information

Code	Description
HealFix	
800500600	HealFix Ti Suture Anchor, 2.0mm, w/two #2-0 sutures W&WB
800500800	HealFix Ti Suture Anchor, 2.0mm, w/two #2-0 sutures W&WB, w/needles
800501400	HealFix Ti Suture Anchor, 2.8mm, w/#2 suture WB
800501200	HealFix Ti Suture Anchor, 2.8mm, w/two #2-0 sutures W&WB
800501000	HealFix Ti Suture Anchor, 2.8mm, w/two #2-0 sutures W&WB, w/needles
800503500	HealFix Ti Suture Anchor, 3.5mm, w/two #2 sutures W&WB
800503300	HealFix Ti Suture Anchor, 3.5mm, w/two #2 sutures W&WB, w/needles
800503900	HealFix Ti Suture Anchor, 5.0mm, w/two #2 sutures W&WB
800503700	HealFix Ti Suture Anchor, 5.0mm, w/two #2 sutures W&WB, w/needles
801310000	1.6mm Drill for 2.0mm HealFix
801310100	Drill Gudie for 2.0mm HealFix
801310200	1.8mm Drill for 2.8mm HealFix
801310300	Drill Guide for 2.8mm/3.5mm HealFix
801310400	2.5mm Drill for 3.5mm/5.0mm HealFix
801310600	3.5mm Drill for 6.5mm HealFix

Code	Description
HealFix Plus	
800504500	HealFix Plus Ti Suture Anchor, 4.5mm, w/two #2 sutures W&WB
800504700	HealFix Plus Ti Suture Anchor, 5.5mm, w/two #2 sutures W&WB
801310800	3.5mm Drill for 4.5mm HealFix Plus
801310900	4.5mm Drill for 5.5mm HealFix Plus
801311000	5.5mm Drill for 6.5mm HealFix Plus
HealFix Mini	
800505300	HealFix Mini Ti Suture Anchor, 1.3mm, w/3-0# suture W, w/tapercut needles
800505400	HealFix Mini Ti Suture Anchor, 1.3mm, w/3-0# suture WB, w/tapercut needles
800505900	HealFix Mini Ti Suture Anchor, 1.8mm, w/2-0# suture W, w/tapercut needles
800506000	HealFix Mini Ti Suture Anchor, 1.8mm, w/2-0# suture WB, w/tapercut needles
801911100	Drill for 1.3mm HealFix Mini
801911200	2.0mm Drill for 1.8mm HealFix Mini